

# Universal Analysis Software v2.0 and above

Reference Guide

VEROGEN PROPRIETARY
Document # VD2019002 Rev. D
May 2021

# **Revision History**

Document#	Date	Description of Change
VD2019002 Rev. D	May 2021	Updated software descriptions to Universal Analysis Software v2.3, which supports the ForenSeq MainstAY Kit with NA24385 as the positive control and refreshes the user interface.  • Added theForenSeq MainstAY assay, Analysis Method, and STR reports.  • Added run metrics for ForenSeq MainstAY libraries.  • Added instructions for performing a status check for troubleshooting.  • Updated instructions for creating CODIS reports, including the addition of options for batch IDs, case IDs, kit names, and partial profiles.  Updated the safety data sheet (SDS) link to verogen.com/product-documentation.  Renamed the software to Universal Analysis Software and this guide to Universal Analysis Software v2.0 Reference Guide.  Added the term stutter to the glossary.  Corrected target run metrics for ForenSeq Kintelligence libraries. Replaced kinship with Kintelligence.
VD2019002 Rev. C	February 2021	Updated software descriptions to ForenSeq Universal Analysis Software v2.2, which supports the ForenSeq Kintelligence Kit with NA24385 as the positive control.  • Added the ForenSeq Kintelligence assay, Verogen Kintelligence Analysis Method, and Kintelligence reports.  • Added run metrics for ForenSeq Kintelligence libraries.  • Added instructions for assigning Unique Dual Index adapters, accessing the analysis folder, switching analysis methods, and printing.  Updated the ForenSeq trademark to a registered trademark.  Refreshed the format of the guide:  • Updated the layout and formatting.  • Updated fonts, table styles, and other design elements.  Added guidance for failed run metrics.  Added a link to the Verogen Documentation web page.  Added the glossary terms clusters, intensity, intralocus balance, positions called, reference SNP ID number, short tandem repeat, typed and untyped, and Unique Dual Index.  Added ambient temperature and reagent chiller temperature as potential sources of high phasing or prephasing values.  Reorganized information on primary workspaces, reports, analysis methods, and reviewing sample data.  Renamed the Sample Details page to Results and updated other user interface terminology as needed.  Clarified that Dark Mode is the only software setting available to non-administrators.  Corrected flow cell compatibility for the ForenSeq mtDNA Control Region Kit.

Document #	Date	Description of Change
VD2019002 Rev. B	August 2020	Updated software descriptions for ForenSeq Universal Analysis Software v2.1, which provides the following enhancements:  • Support for the ForenSeq mtDNA Whole Genome Kit  • A Run, Project, and Sample History setting in System Settings  • Project History and Sample History worksheets for Project and Sample Reports Updated guide format and organization, including adding an index and consolidating workspace descriptions.  Revised step-by-step instructions to be more succinct and include field descriptions. Updated telephone numbers for Verogen contact information.  Updated glossary term definitions and added the terms coverage plot, indel, library, mtDNA navigator, plexity, position viewer, system event, and user action.  Updated figures to include call-outs highlighting user interface elements.  Added information on strand bias and lost passwords.  Added instructions for creating projects, assigning a sample to multiple projects, modifying user-created analysis methods, and enabling and disabling dark mode.  Added another possible cause of low cluster density: belatedly combining bead-based normalized libraries and Hybridization Buffer (HT1).  Noted that the analysis folder contains VCF and BAM files.  Noted that the software follows SWGDAM nomenclature rules.  Noted the conditions for displaying the Call or No Call filter.  Clarified that each sample name must be unique in the software.  Clarified that phasing and prephasing values are for Read 1 only.
VD2019002 Rev. B	August 2020	Distributed troubleshooting information among individual sections.  Moved information on analysis method settings to <i>Analysis Defaults for ForenSeq Universal Analysis Software v2 (Pub. No. VD20200045)</i> .  Corrected duplicate entries of D7S820 and D16S539 in the Human Sequencing Control. Corrected the list of information included in project and sample reports.  Corrected citations that appear in the appendix.
VD2019002 Rev. A	August 2019	Initial release

# **Table of Contents**

Revision History	
Table of Contents	4
Chapter 1 Software Overview	7
Introduction	7
Accessing the Software	
Primary Workspaces	
Software Status	
Glossary	
Chapter 2 Run Setup	25
Create a Run	
Edit a Run	
Add Samples to a Run	
Manage Samples in a Run	
Monitor a Run	
Chapter 3 Data Analysis	30
Review Run Metrics	30
Review mtDNA Results	31
Review Kintelligence Results	33
Review MainstAY Results	34
Modify Kintelligence and MainstAY Results	35
Add Analyzed Samples to a Project	
Reanalyze Samples	37
Chapter 4 Report Generation	39
Generating mtDNA Reports	39
Generating Kintelligence Reports	
Generating MainstAY Reports	
Chapter 5 Run, Project, and Sample Management	49
Create a Project	AC

Edit a Project	49
Obtain the Analysis Folder Path	50
Printing	50
History	51
Chapter 6 Software Settings	55
Settings Overview	55
Analysis Methods	
CODIS Report Defaults	
User Management	
System Settings	61
Enable or Disable Dark Mode	62
Chapter 7 Troubleshooting	63
Troubleshoot a Run	
Troubleshoot Analysis Errors	
Perform a Status Check	
- Griorin a Grada Griodic	
Appendix A Run Metrics	68
Accessing Run Metrics	68
Quality Metrics	68
Positive and Negative Controls	73
Sample Representation	75
Appendix B Report Formats	78
File Formats for mtDNA Reports	78
File Formats for Kintelligence Reports	
File Formats for MainstAY Reports	
Appendix C Supporting Information	85
Special Handling of mtDNA Data	
Special Handling of MiDNA Data	
Gender Estimation for Kintelligence Samples	
Contributor Status for MainstAY Samples	
Gender Estimation for MainstAY Samples	

Stutter filtering for MainstAY Samples	. 90
Index	91
Technical Support	.95

This document and its contents are proprietary to Verogen, Inc. and its affiliates, and are intended solely for the contractual use of its customer in connection with the use of the product(s) described herein and for no other purpose. This document and its contents shall not be used or distributed for any other purpose and/or otherwise communicated, disclosed, or reproduced in any way whatsoever without the prior written consent of Verogen. Verogen does not convey any license under its patent, trademark, copyright, or common-law rights nor similar rights of any third parties by this document.

The instructions in this document must be strictly and explicitly followed by qualified and properly trained personnel in order to ensure the proper and safe use of the product(s) described herein. All of the contents of this document must be fully read and understood prior to using such product(s).

FAILURE TO COMPLETELY READ AND EXPLICITLY FOLLOW ALL OF THE INSTRUCTIONS CONTAINED HEREIN MAY RESULT IN DAMAGE TO THE PRODUCT(S), INJURY TO PERSONS, INCLUDING TO USERS OR OTHERS, AND DAMAGE TO OTHER PROPERTY.

VEROGEN DOES NOT ASSUME ANY LIABILITY ARISING OUT OF THE IMPROPER USE OF THE PRODUCT(S) DESCRIBED HEREIN (INCLUDING PARTS THEREOF OR SOFTWARE).

© 2021 Verogen, Inc. All rights reserved.

All other names, logos, and other trademarks are the property of their respective owners.

For Research, Forensic, or Paternity Use Only. Not for any clinical or therapeutic use in humans or animals.

# Software Overview

ntroduction	7
Accessing the Software	8
Primary Workspaces	9
Software Status	19
Glossary	20

#### Introduction

Universal Analysis Software (UAS) v2.0 is a forensic genomics solution for run setup, sample management, data analysis, and reporting. Preinstalled on a dedicated server, the software integrates with the MiSeq FGx<sup>®</sup> Sequencing System to automatically generate FASTQ files, perform alignment, and call variants from run data.

The software analyzes three types of sequencing data for the following workflows:

- mtDNA analysis—For libraries prepared with the ForenSeq® mtDNA Control Region Kit, the software interrogates
  the control region of the mitochondrial genome (mtGenome). For ForenSeq mtDNA Whole Genome libraries, the
  software interrogates the entire mtGenome. UAS v2.0 also analyzes mitochondrial DNA (mtDNA) data from
  libraries prepared with custom or third-party assays and sequenced on the MiSeq FGx System.
- **Kintelligence analysis**—For libraries prepared with the ForenSeq Kintelligence Kit, the software interrogates single nucleotide polymorphisms (SNPs) for upload to GEDmatch<sup>®</sup> PRO.
- **MainstAY analysis**—For libraries prepared with the ForenSeq MainstAY Kit, the software interrogates autosomal short tandem repeats (auSTRs) and Y-STRs.

The analysis method varies for each library type and the reports are different. However, run setup is the same.

# Compatibility of Run Components

Use the version of UAS v2.0 that is compatible with your library prep kit and a flow cell that supports the desired plexity and read coverage. The MiSeq FGx Reagent Kit provides the standard flow cell and the MiSeq FGx Reagent Micro Kit provides the micro flow cell.

Table 1 Version compatibility

Library Prep Kit	Flow Cell	Maximum Plexity*	Software Version
ForenSeq MainstAY Kit	Micro	96 libraries	ForenSeq UAS v2.3, or later
ForenSeq Kintelligence Kit	Standard	3 libraries	ForenSeq UAS v2.2, or later
ForenSeq mtDNA Whole Genome Kit	Standard	16 libraries	ForenSeq UAS v2.1, or later
ForenSeq mtDNA Control Region Kit	Micro	48 libraries	ForenSeq UAS v2.0, or later

<sup>\*</sup> When using the default analysis method

#### **Analysis Workflow**

The following diagram outlines the analysis workflow, which includes UAS v2.0 and two applications that run on the MiSeq FGx System: MiSeq FGx Control Software (MFCS) and Real-Time Analysis (RTA). For more information on these applications, see the MiSeq FGx Sequencing System Reference Guide (document # VD2018006).

#### Figure 1 Overview of data analysis steps



#### Sign In or Create Account

The user enters a user name and password to access the software.



#### Set Up a Run

The user manually inputs sample information or imports a .txt file.



#### Sequence Libraries

The control software images the flow cell and RTA performs image analysis and base calling.



#### **Analyze Data**

The software uses data from RTA to align reads and make variant calls with quality indicators. The user reviews the results in the software.



#### **Generate Reports**

The user generates reports on the analysis results in a variety of formats.



UAS v2.0



MiSeq FGx System

# Accessing the Software

UAS v2.0 is a browser-based application that does not require an internet connection. Open a web browser (Chrome is recommended) on a computer with access to the same network as the MiSeq FGx System. Set screen resolution to at least 992 pixels for optimal display.

To access the UAS v2.0 user interface, enter the server address in the web browser address bar. If the Domain Name System (DNS) does not list the server, enter the server IP address for the network.

#### Create an Account

Before using UAS v2.0, you must create an account. The first user to create an account is automatically assigned administrator privileges. The administrator is responsible for approving new accounts. For instructions, see *Enable and Disable Accounts* on page 61.

Lost passwords are not recoverable, so creating a second administrator account is highly recommended. If the administrator password is lost and only one administrator account exists, the software must be uninstalled and reinstalled.

- 1. Access the network server through a web browser.
- 2. On the Login page, select Register as a new user.
- 3. On the Register page, complete the following fields:
  - Email—Enter an email address to be the user name for the account.
  - Password

    —Enter a password for the account.
  - Confirm password—Reenter the password for the account.
- 4. Select **Register** to submit your account for approval from the system administrator.

A confirmation message appears. You can access the software after approval.

#### Sign In or Sign Out

- 1. Sign in as follows.
  - a. Access the network server through a web browser.
  - b. Enter your email address (user name) and password.
  - c. Select Login.

The software opens to the Home page.

2. When you are ready to sign out, select Account, and then select Logout.

# Change a Password

- 1. On the navigation rail, select **Account**.
- 2. Select Change Password.
- 3. In the Current Password field, enter your password.
- 4. In the New Password field, enter a new password.
- 5. In the Confirm Password field, repeat the new password.
- 6. Select **Save** to apply the new password.

# **Primary Workspaces**

The following pages comprise the primary workspaces in UAS v2.0:

Home—View and create runs and projects. This page functions as a system dashboard.

- Runs—Set up and manage runs before and during a run.
- Projects—View and manage projects and samples after a run and generate reports.
- Results—Review data on a sample-by-sample basis after a run.


The software also includes printing, commenting, and other functions that supplement run setup, data analysis, and reporting. For instructions on using these supplementary features, see *Run, Project, and Sample Management* on page 49.

#### Home Page and Navigation

The Home page displays information and functions organized into runs and projects. Selecting a run goes to the Run Details page and selecting a project goes to Projects. Color-coded Quality Metrics icons provide an overall view of run performance. Hovering over an icon displays more detailed metrics.

The navigation rail vertically arranges icons and puts key destinations within easy reach. It appears consistently on the left side of each page so you can go to runs, projects, settings, and your account or return to the Home page.

Figure 2 Overview of the Home page

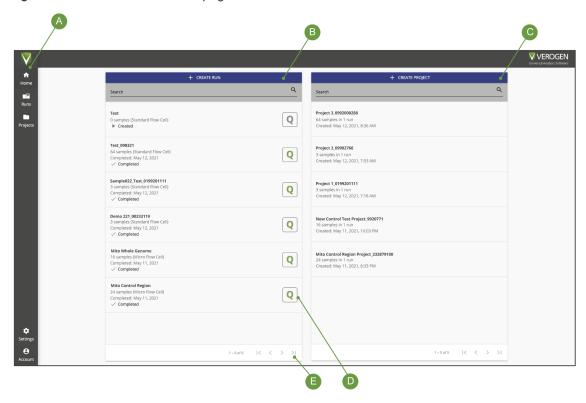

- A Navigation rail—Icons that go to key destinations
- B Runs-List of runs with quality metrics and options to search and create
- C **Projects**—List of projects with options to search and create
- D Quality Metrics icon—Quality metrics for a completed run
- E Scroll arrows-Move through all runs or projects

#### Runs


Runs is the sequencing workspace for creating and monitoring runs and the input and organization of sample information. Selecting Runs on the navigation rail opens the Runs navigation drawer, which provides functions for creating and searching runs.

The Runs navigation drawer also lists every run in the software with the run state, number of samples in the run, flow cell type, and color-coded Quality Metrics icon. Runs are listed from newest to oldest. Use the scroll arrows to move through the complete list.

#### Run Details Page

Selecting a run from the Runs navigation drawer opens the Run Details page. Each Run Details page is a run-specific workspace for managing the run and samples in the run. Throughout the software, selecting a run name hyperlink goes to the Run Details page.

Figure 3 Run Details Page on the Universal Analysis Software user interface

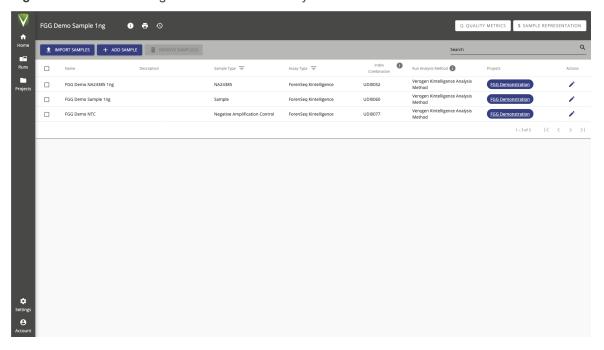

A table displays all samples in the run with information about each sample and options for sorting and filtering. Functions for importing, adding, searching, and deleting samples appear above the table. The header section of a Run Details page displays the run name, analysis method, quality metrics, sample representation, and the following icons.

Table 2 Icons on a Run Details page

Icon	Name	Function
0	Run Details	Edit the run and view the run state, cycles per read, and other run information.
<b>3</b>	Run Activity	See system events and user actions for the run and add comments.  This icon is visible when the history setting is enabled. See <i>System Settings</i> on page 61.
6	Print	Print the Run Details page.

# **Projects**

Projects is the post-sequencing workspace for creating and managing projects, accessing analysis results, and generating reports. Selecting Projects on the navigation rail opens the Projects navigation drawer, which provides functions for creating and searching projects.

The Projects navigation drawer also lists every project in the software with the number of samples and runs in the project and the date and time of project creation. Projects are listed from newest to oldest. Use the scroll arrows to move through the complete list.

#### **Project Side Sheet**

Selecting a project from the Projects navigation drawer opens the Project side sheet and the Samples page. The Project side sheet lists all samples in the project, grouped by run. Expand the runs to see which samples are included. A green checkmark indicates an analyzed sample.

The Project side sheet provides the following icons for accessing additional pages and project information.

Table 3 Icons on a Project side sheet

Icon	Name	Function
1	Edit	Change the project name and add or change the optional description.
8	Samples	Open the Samples page for the project.
	Reports	Open the Reports page for the project.
<b>1</b>	Project Activity	See system events and user actions for the project and add comments.  This icon is visible when the history setting is enabled. See <i>System Settings</i> on page 61.
ō	Print	Print the Samples page.
QPN	QPN	Open the QPN page for each run in the project. Each icon letter is color-coded to provide high-level status. See <i>Run Metrics</i> on page 68.

#### Samples Page

The Samples page lists information for each sample in a selected project with functions for adding, reanalyzing, and searching samples. Filters narrow results by sample or assay type and sort by sample name, run name, and analysis method.

A QC Indicators column displays the number of quality control (QC) indicators for a sample. N/A indicates a sample with zero QC indicators. Samples with QC indicators require further review.

#### Reports Page

The Reports page lists reports in a project by creation date, starting with the most recent report. On this page, you can create reports or find and download existing reports. Sort the Date Created and File Name columns and filter the Report Type column to organize results.

#### **QPN Page**

The QPN page is divided into three tabs: Quality Metrics, Positive and Negative Controls, and Sample Representation. Together, these tabs provide an overall view of run performance. For details, see *Run Metrics* on page 68.

#### Results Page for mtDNA

Selecting an mtDNA sample from a project opens the Results page. Each Results page is a sample-specific workspace for in-depth analysis of results. A header section displays the sample name, analysis method applied to the sample, buttons to create reports and filter and compare samples, and the following icons.

Table 4 Icons on an mtDNA Results page

Icon	Name	Function
6	Sample Details	View the sample type and analysis state, analysis folder path, and other analysis information.
<b>3</b>	Sample Activity	See system events and user actions for the sample and add comments.  This icon is visible when the history setting is enabled. See <i>System Settings</i> on page 61.

The remainder of the Results page is divided into an mtDNA navigator, position viewer, and coverage plot that examine each position. Changing the position in one of these tools simultaneously updates the position in the others.

#### mtDNA Navigator

The mtDNA navigator displays data spanning the whole mtGenome or control region. Data include the IT, variants, calls or no calls, and coverage at various positions. The movable pointer updates the position.

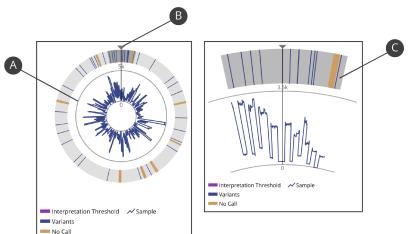


Figure 4 mtDNA navigators for whole genome (left) and control region (right) results

- A Maximum coverage indicator
- B Movable pointer
- C Control region view (filter applied)

By default, the mtDNA navigator shows the entire mtGenome. To focus on the control region only, apply the Control Region filter.


#### **Position Viewer**

The position viewer compares the sample sequence to the revised Cambridge Reference Sequence (rCRS). Color-coded dashes indicate ITs, variants, and no calls. For example, the blue dash in Figure 5 underlines the base at position 16,362 because the call (C) differs from the reference (T).

Individually select or enter each position number, or use the scroll arrows to move through pages of sequence segments and hotspots. For each selected position, the position viewer displays the following data:

- Reference and call—The base calls for the reference and sample.
- Total read count—The total number of reads from the forward and reverse strands.
- Strand depth—The strand (forward or reverse) read count with the majority of reads.
- Percentages—The percentage of each base compared to the AT and IT.

Figure 5 Position viewer set to position 16,362

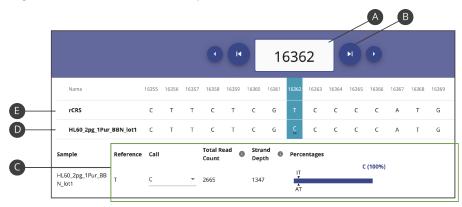

- A Selected position number
- B Scroll arrows
- C Data for the selected position
- D Sample sequence
- E Reference sequence


#### Strand Bias

If you want to evaluate strand bias, which provides an additional metric for evaluating quality, use the strand depth and total read count. The software reports strand depth as the read count from the forward or reverse strand that has the majority of reads for a base position. Subtract from the total read to determine the read count for the opposite strand.

#### Coverage Plot

The coverage plot displays the coverage (number of reads) for each mtDNA position. Positions are plotted along the x-axis and the y-axis shows coverage. A pointer marks the selected position. Move the pointer to shift the position.

Figure 6 Coverage at each position

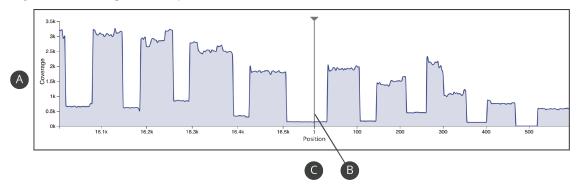

- A Coverage along the y-axis
- B Movable pointer
- C Positions along the x-axis

#### Results Page for Kintelligence

Selecting a Kintelligence sample from a project opens the Results page. Each Results page is a sample-specific workspace for in-depth analysis of results. A header section displays the sample name, analysis method applied to the sample, button to create reports, and the following icons.

 Table 5
 Icons on a Kintelligence Results page

Icon	Name	Function
0	Sample Details	View the sample type and analysis state, analysis folder path, and other analysis information.
<b>3</b>	Sample Activity	See system events and user actions for the sample and add comments.  This icon is visible when the history setting is enabled. See <i>System Settings</i> on page 61.
<b>6</b>	Print	Print the Sample Overview section of the Results page.

In addition to the header, the Results page for a Kintelligence sample includes a Sample Overview section and a Sample Details section.

#### Sample Overview Section

The Sample Overview section displays locus call rate, total reads, gender, contributor status, and heterozygote and intralocus balance metrics. Under these data, a table displays the percentage and total number of typed loci relative to base pair length per SNP type. The entire section is collapsible and you can edit both gender and contributor status.

#### Sample Details Section

The Sample Details section displays data for the SNP types listed in the following table. A card represents each locus, displaying the locus name, chromosome, allele, intensity, and intralocus balance. Each allele is color-coded: A is green, T is red, C is blue, and G is gray.

SNP Type	Description
X-SNP	Informs lineage and estimates gender in combination with Y-SNPs.
Y-SNP	Informs lineage and estimates gender in combination with X-SNPs.
Ancestry SNP (aSNP)	Estimates biogeographical ancestry.
Phenotype SNP (pSNP)	Estimates hair and eye color.
Identity SNP (iSNP)	Informs on source attribution.
Kinship SNP (kSNP)	Informs on potential biological associations (relatives).

#### Locus Details

Selecting a card in the Sample Details section opens locus details, which include allele intensity, thresholds, intralocus balance, and quality control (QC) indicators. Under Typed, a toggle key indicates whether the allele is typed.

- The Previous Locus and Next Locus arrows move through details for each locus.
- The Collapse icon closes locus details.

#### Locus Quality Indicators

Each card displays a color-coded horizontal bar: gray indicates a locus with no QC indicators and orange indicates a locus with at least one QC indicator. Any applicable indicators are displayed on the card.

Table 6 QC indicators

Symbol	Name	Description
	Allele Count	The locus has more typed alleles than expected for a single-source sample of the indicated gender.
?	Unexpected Allele	The locus has an unexpected allele (non-reference or alternate) that is greater than or equal to the AT.
90	Imbalanced	The locus is not homozygous, not hemizygous, and falls below the intralocus balance setting.
O	Analytical Threshold	The locus has expected allele signal below the AT and no alleles greater than or equal to the IT.
0	Interpretation Threshold	The locus has expected allele signal greater than or equal to the AT, but below IT.
$\otimes$	Not Detected	The software did not detect signal for the locus.
	User Modified	Someone edited the typed status of an allele at the locus.

# Results Page for MainstAY

Selecting a MainstAY sample from a project opens the Results page. Each Results page is a sample-specific workspace for in-depth analysis of results. A header section displays the sample name, analysis method applied to the sample, button to create reports, and the following icons.

Table 7 Icons on the MainstAY Results page

Icon	Name	Function
0	Sample Details	View the sample type and analysis state, analysis folder path, and other analysis information.
<b>5</b>	Sample Activity	See system events and user actions for the sample and add comments.  This icon is visible when the history setting is enabled. See <i>System Settings</i> on page 61.

In addition to the header, the Results page for a MainstAY sample includes a Sample Overview section and detailed coverage and histogram views.

#### Sample Overview Section

The Sample Overview section displays locus call rate, total reads, gender, and contributor status. You can modify both the gender and contributor status of the sample.

#### Coverage View

The coverage view displays coverage by loci or alleles. By default, loci coverage is displayed and depicted as a heatmap. Divided into auSTRs and Y-STRs, the heatmap shows the number of loci typed out of the total number of possible loci. Hovering over a green segment displays the locus name, amplicon length in base pairs (bp), and intensity. Gray segments do not contain these data.

Switching the view to alleles opens a scatter diagram that plots amplicon length along the X-axis and intensity along the Y-axis. Each dot represents a typed allele in the sample: yellow dots indicate auSTRs and blue dots indicate Y-STRs. Hovering over a dot reveals the STR type, locus name and allele number, amplicon length, and intensity value. Selecting the portion of the allele scatterplot using a screen-select allows a user to zoom into this section of the visualization. Clicking on a selected locus square on the heatmap or on an individual dot on the allele scatter diagram takes you directly to that locus in the histogram.

#### Histogram View

The histogram view displays intensity values and thresholds for each locus in the sample. All loci are displayed by default, but you can switch to viewing only auSTR loci or only Y-STR loci. You can also use the filtering, sort, and search options to refine results.

Coverage plots can be hidden by clicking the **Hide Coverage Plots** button on the top right corner of the user interface. Hiding the coverage plot scales the histogram view to the larger screen area. Collapsing the filter bar or the project navigation drawer by selecting the grey bar with the arrow also scales up the histogram view. Each locus is depicted as a card that displays a bar chart with one bar for each allele. Hovering over a bar reveals the intensity. Effective analytical thresholds appear as gray shading so you can confirm whether each allele meets the thresholds.

#### Locus Details

Selecting a card opens locus details, which enlarges the bar chart and shows additional data for each allele, including intralocus balance, thresholds, intensity, stutter, sequences, and QC indicators. Toggle keys indicate whether alleles are typed.

- The Previous Locus and Next Locus arrows move through details for each locus.
- The X icon closes locus details.

#### Locus Quality Indicators

Each card displays a color-coded horizontal bar: green indicates a locus with no QC indicators and orange indicates a locus with at least one QC indicator.

Table 8 QC indicators

Symbol	Name	Description
	Allele Count	The locus has more typed alleles than expected for a single-source sample of the indicated gender.
9	Imbalanced	The locus is not homozygous, not hemizygous, and falls below the intralocus balance setting.
O	Analytical Threshold	The locus has expected allele signal below the AT and no alleles greater than or equal to the IT.
O	Interpretation Threshold	The locus has expected allele signal greater than or equal to the AT, but below IT.
$\otimes$	Not Detected	The software did not detect signal for the locus.
	User Modified	Someone edited the typed status of an allele at the locus.
	Stutter	At least one sequence exceeds the stutter filter percentage at a stutter position of a possible parent allele.

# Software Status

The software displays a state for runs and analyses. The state determines which modifications are permitted and are also useful for general status updates.

#### Run State

After a run is created and saved, you can view the run state by selecting the Information icon on the Run Details page displays the run state. The following table describes each possible state.

State	Description
Created	The run information is saved in UAS v2.0.
Sequencing	The run is in process on the MiSeq FGx System.
Completed	The run completed successfully.
Errored	An error prevented the run from finishing or it was manually stopped.
Paused	A user temporarily suspended the run.
Stopped	A user ended the run due to incorrect setup, poor data quality, or other issue.

# **Analysis State**

After a sample is added to a run, you can view the analysis state by selecting the Information icon on the Results page. The following table describes each possible state.

State	Description
Created	UAS v2.0 has not started analysis.
Queued	The analysis is in line to start.
Processing	Analysis or reanalysis is in progress.
Completed	Analysis has successfully completed.
Errored	Analysis has unsuccessfully completed.

# Glossary

Term	Description
Amplicon	The product of PCR amplification of a targeted region of interest in a genome.
Amplicon end position	The last base position of an amplicon, including the PCR primer sequence.
Amplicon start position	The first base position of an amplicon, including the PCR primer sequence.
Analysis method	Settings and thresholds that inform the analysis of sequencing data to generate allele or variant calls and QC indicators.
Analytical threshold (AT)	The percentage defined in an analysis method that a read count must reach for the software to call a base or type an allele. Signal below the AT might be visible, but is not typed.
Assay type	The library prep kit used to process samples for sequencing and analysis.
Call	The reported base for a position (coordinate) based on the analysis method.
Cluster density	The number of clusters on a flow cell per square millimeter (K/mm²).
Clusters	A clonal group of DNA bound to the surface of a flow cell that produces one read. One template DNA strand seeds each cluster and is clonally amplified using bridge amplification until the cluster has about 1000 copies.
Clusters passing filter	The percentage of clusters that pass quality filters. For more information, see <i>Clusters Passing Filter</i> on page 69.
Control region	A subset of the human mtGenome (positions 16,024-576) containing three hypervariable regions.
Coverage	The average number of reads for a locus, allele, or variant position in a sample that align to or cover known reference bases. Sequencing coverage level (read depth) can help determine whether an allele or variant call is made with confidence.
Coverage plot	A graph on a Results page that displays the coverage at each position in an mtDNA sample. For details, see <i>Coverage Plot</i> on page 15.
Cycle	A chemistry step and an imaging step to call one base in a sequencing run. For cycle numbers, see <i>Cycles per Read</i> on page 70.

Term	Description
Deletion	A base deleted from the sample sequence that is present in the reference sequence.
Differences indicator	A software feature that highlights different allele or variants calls between two samples being compared.
Flow cell	A glass slide coated with oligos complementary to the library adapter sequences, allowing the library to adhere for a run.
Forward primer length	The number of bases in the forward PCR primer (oligonucleotide).
Forward primer start coordinate	The starting base at the $5^{\prime}$ end of the forward PCR primer. The base depends on the amplicon and analysis method.
Hotspot	An mtDNA position with a variant, interpretation threshold, differences, no call, call indicator, or user modification.
Human Sequencing Control (HSC)	A pre-amplified and tagged positive control added to a library pool to facilitate troubleshooting of instrument and sequencing reagent issues.
Index adapter sequence	A short stretch of oligonucleotides added to each sample during library prep for identification purposes. The software uses the oligonucleotide sequences to demultiplex pooled libraries.
Index CV	The percentage of reads assigned to each library, representing the distribution of read counts of the libraries in the run. CV is the coefficient of variation for the number of read counts across all indexes.
Indel	The insertion or deletion of bases in a genome.
Insertion	A base in the sample sequence that is not present in the reference sequence.
Intensity	The software reports signal intensity as the number of reads.
Intralocus balance	The balance of read counts between typed alleles at a heterozygous locus. The balance is measured as the intensity of the minimum intensity typed allele divided by the intensity of the maximum intensity typed allele. When intensities are identical, the intralocus balance is 100%.
Interpretation threshold (IT)	The percentage defined in an analysis method that determines the read count that must be reached for the software to call a base or type an allele. Signal below the IT but greater than or equal to the AT can be manually typed.
Library	A DNA sample that uses primers and adapters to ligate specific regions of the genome for sequencing.
Library type	The function of a library, either primer-directed sequencing or non-directed sequencing, for custom assays and analysis.
Minimum quality score	The Q-score reads must meet to contribute to the variant call at each DNA coordinate. This value is defined in an analysis method.
Minimum read count	The minimum read count a position must reach for a base to be called.

Term	Description
mtDNA navigator	A section of a Results page that displays variants and other sequencing results for an mtDNA sample. For details, see <i>mtDNA Navigator</i> on page 13.
N indicator	The quality indicator for reagent blanks and negative amplification controls.
Negative amplification control	A sample that does not contain DNA as a template for amplification.
No call	The call assigned when a position lacks read counts for a specific base.
Original call	An unmodified variant call based on the run analysis method.
P indicator	The quality indicator for positive amplification controls that use control DNA, such as HL60.
Percent A, C, G, T	The percentage of forward and reverse reads out of the total reads for the A, C, G, or T bases at the call position.
Percent del	The percentage of forward and reverse reads out of the total reads for a deletion.
Percent ref	The percentage of forward and reverse reads out of the total reads that do not support an insertion at an insertion point in mtDNA.
Percentages	A position viewer feature that shows the percentage of each nucleotide called at a position in relation to the thresholds.
Phasing	The percentage of bases in a cluster that fall behind the current cycle in Read 1. For details, see <i>Phasing and Prephasing</i> on page 69.
Plexity	The number of libraries combined in a pool for sequencing. For example, if eight libraries are combined in a pool, the plexity is eight.
Position	A base coordinate of mtDNA using the numbering convention of the rCRS, comprised of positions 1-16,569.
Position viewer	A section of a Results page that displays reference and sample sequences and data for the selected position in an mtDNA sample. For details, see <i>Position Viewer</i> on page 14.
Positions called	How many positions out of the total number of interrogated mtDNA positions, including indels, the software called a base for.
Positive amplification control	A known sequence of DNA provided in the library prep kit and prepared with forensic and reference samples to ensure that library prep and sequencing reagents are functional.
Prephasing	The percentage of bases in a cluster that jump ahead of the current cycle in Read 1. For details, see <i>Phasing and Prephasing</i> on page 69.
Project	A collection of analyzed results for at least one sample, and the primary workspace for viewing and modifying results and creating reports.
Revised Cambridge Reference Sequence (rCRS)	The reference genome for calling variants in the mtGenome.
Read	The sequence string of A, T, C, and G bases corresponding to an amplicon. A run can simultaneously generate millions of reads from multiple loci and multiple samples.

Term	Description
Reads A, C, G, T	The number of forward and reverse reads for the A, C, G, or T base at the DNA call position.
Reads del	The number of forward and reverse reads for a deletion at the call position.
Reads ref	The number of forward and reverse reads that do not support an insertion at an insertion point.
Reagent blank	An extraction control without DNA that is carried through the entire sample processing.
Read 1 (i5) primer sequence	The sequence of the PCR primer that is tagged with the Index 2 (i5) adapter sequence. Index 2 is part of the adapter sequence, as is the binding site for the Read 1 sequencing primer.
Read 2 (i7) primer sequence	The sequence of the PCR primer that is tagged with the Index 1 (i7) adapter sequence. Index 1 is part of the adapter sequence, as is the binding site for the Read 2 sequencing primer.
Reagent cartridge kit lot	The batch code to identify the batch or lot that a component of the MiSeq FGx Reagent Kit or MiSeq FGx Reagent Micro Kit was manufactured in.
Reanalysis	The process of reanalyzing a sample with another analysis method, after the original run analysis method.
Reference call	The base call of the rCRS at a given position.
Reference SNP ID number (rsID)	An identification tag the National Center for Biotechnology Information (NCBI) assigns to a group of SNPs that map to the same location.
Reverse primer length	The number of bases in the reverse PCR primer (oligonucleotide).
Reverse primer start coordinate	The starting base at the 5' end of the reverse PCR primer. The base depends on the amplicon and analysis method.
Run	The process of clustering and sequencing libraries to generate base calls.
Run analysis method	The analysis method assigned to a sample for the initial run analysis.
Sample	The DNA input for a library prep kit, before the addition of primers and adapters.
Sample compare	A software feature that compares an mtDNA sample of interest to other mtDNA samples.
Sample of interest	A user-designated sample that is compared to other selected samples.
Sample plate	The 96-well plate that contains libraries for sequencing.
Sample representation	The number of reads per sample for a run. These data are available after sequencing and analysis are complete.
Sample type	The function of a sample: forensic or reference sample, positive amplification control, negative amplification control, and reagent blank.
Sample well position	The column and row position in a 96-well plate containing samples or libraries. For example, position A1 indicates a well at the intersection of row A and column 1.

Term	Description
Short tandem repeat (STR)	Areas of the genome where short sequences of DNA (2-6 bp) are repeated, such as (GATA)n.
Single-nucleotide polymorphism (SNP)	Variation of one base pair in a DNA sequence relative to a known DNA reference sequence.
Strand depth	The read count from the forward or reverse strand that has the majority of reads for a base position.
Strand direction	The selection of a forward or reverse strand for a primer-directed library type.
Stutter	Polymerase slippage that can occur during PCR amplification of repetitive DNA sequences, library prep, cluster generation, or sequencing. The slippage can create DNA amplification products that are less than or greater than the size of the parent allele.
System event	An action the software automatically performs and records.
Total read count	The total number of reads from the forward and reverse DNA strands for a base position.
Typed	An allele call that contributes to the locus genotype.
Unique Dual Index (UDI)	Premixed Index 1 (i7) and Index 2 (i5) adapters for dual index combinations that have no redundant sequences.
Untyped	A no-call allele that does not contribute to the locus genotype.
User action	A user interaction with the software, such as adding a sample to a run.
User-modified call	A base call that a user changed.
Variant	A base at a position in an mtDNA sequence that differs from the rCRS base and a synonym for SNP.
Whole mitochondrial genome	The entire 16,569 bp mtGenome.

# Run Setup

Create a Run	25
Edit a Run	25
Add Samples to a Run	25
Manage Samples in a Run	27
Monitor a Run	29

#### Create a Run

- 1. On the Home page, select Create Run.
- 2. In the Create Run dialog box, complete the following fields:
  - Run Name

    Enter a preferred name to identify the run.
  - Description—Enter a description of the run.
  - Flow Cell Type—Select Standard or Micro, depending on your reagent kit.
  - Reagent Cartridge Kit Lot #-Enter the LOT number printed on the reagent cartridge packaging.

The user interface displays required fields with an asterisk (\*).

3. Select Save to create the run.

The Run Details page appears.

#### Edit a Run

- 1. On the navigation rail, select **Runs**, and then select a run.
- Select Information, and then select Edit.
- 3. In the Edit Run dialog box:
  - Edit the run name.
  - Add or modify the optional run description.
  - Select a new flow cell type, Standard or Micro.
  - Add or modify the optional reagent cartridge kit lot number.

Editing the flow cell type requires a run state of Created.

4. Select Save.

# Add Samples to a Run

Add samples and sample information to a run using one of two methods: import a sample sheet or individually add samples. After samples are added, the Run Details page lists all samples included in the run.

All samples in a run must have the same assay type.

When adding samples to a run, you associate each sample with at least one project. If you enter the name of a new project, the software creates the project. Because a sample can be assigned to multiple projects, a run can include samples from multiple projects.

#### Import a Sample Sheet

- 1. On the navigation rail, select **Runs**, and then select a run.
- Select Import Samples.
- [Optional] Select Sample Sheet Template to download a .txt file that provides examples of Sample Sheet formats for each Assay format.
- 4. Select Drop Files Here, and then navigate to the sample sheet location.
- Select the sample sheet (.txt file), and then select **Open**.
   The sample sheet appears under the file upload area. Uploading another .txt file overrides the current file because the software accepts only one sample sheet per run.
- Select Next.
- 7. Review the information from the uploaded sample sheet to ensure accuracy. Scroll to see all entries. If utilizing the template, make sure to delete the sections pertaining to assay types that are not utilized in the current analysis.
- 8. Select **Import** to add the samples to the run.

#### Assigning Unique Dual Indexes

The samplesheet template includes two columns for entering index adapter names: Index 1 (i7 index) and Index 2 (i5 index). When using Unique Dual Index (UDI) adapters, either enter the index name in the Index 1 column or both columns.

#### **Assigning Projects**

The InitialProjectList column of a sample sheet specifies which project a sample is assigned to. To assign multiple projects to a sample, enter each project name separated by a comma followed by a space and put quotes around the entry. For example:

```
"Project 1, Project 2, Project 3"
```

# Individually Add Samples

- 1. On the navigation rail, select **Runs**, and then select a run.
- Select Add Sample.
- 3. In the Add Sample dialog box, complete the following fields:

Field	Instruction
Sample Name	Enter a unique identifier for the sample. Each sample name must be unique in UAS v2.0, regardless of sample type.
Description	Enter a preferred description of the sample.

Field	Instruction
Sample Type	Select a sample type:  • Sample—A forensic or reference sample  • ForenSeq Positive Amplification Control—The control DNA (NA24385, or HL60)  • Negative Amplification Control—The negative control (water)  • Reagent Blank—The extraction control
Assay Type <sup>1</sup>	Select the library prep method:  • ForenSeq mtDNA Whole Genome—Libraries prepared with the ForenSeq mtDNA Whole Genome Kit  • ForenSeq mtDNA Control Region—Libraries prepared with the ForenSeq mtDNA Control Region Kit  • mtDNA Custom—Libraries prepared with a custom or third-party assay  • ForenSeq Kintelligence—Libraries prepared with the ForenSeq Kintelligence Kit  • ForenSeq MainstAY—Libraries prepared with the ForenSeq MainstAY Kit
Index 1 <sup>2</sup>	Select the name of an i7 index adapter.
Index 2 <sup>2</sup>	Select the name of an i5 index adapter.
Index Combination <sup>2</sup>	Select the name of a premixed i7 and i5 index adapter (UDI).
Run Analysis Method	Select a method for analyzing samples in the run.
Initial Projects	Enter the name of a project to assign the sample to. To assign the sample to multiple projects, enter multiple project names.
Sample Plate	Enter a preferred name for the 96-well plate.
Sample Well Position	Enter the position of the plate well the sample occupies.

<sup>&</sup>lt;sup>1</sup> Because the software allows one assay type per run, the assay type assigned to the first sample is automatically assigned to the others.

An asterisk (\*) on the user interface indicates a required field.

4. Select Add to add the sample to the run.

# Manage Samples in a Run

After adding samples to a run, you can edit sample information, change project assignments, or remove samples from the run. Use filters, sorting, and the search function to quickly find samples.

# Find a Sample

- 1. On the navigation rail, select **Runs**, and then select a run.
- 2. To search for samples, enter any of the following information in the Search field:
  - Sample name
  - Description
  - · Index adapter name
  - · Run analysis method

<sup>&</sup>lt;sup>2</sup> Depending on assay type, the software displays either the Index 1 and Index 2 fields or the Index Combination field.

- 3. To filter the samples, do as follows.
  - a. In the Sample Type or Assay Type column, select Filter.
  - b. Select at least one checkbox.
  - c. Reselect Filter to apply the selected filters.
- 4. To sort the samples, select a column heading.

The Name, Description, and Index Combination columns are sortable.

# **Edit Sample Information**

- 1. On the navigation rail, select **Runs**, and then select a run.
- 2. In the Actions column, select Edit Sample Details.
- 3. In the Edit Sample dialog box, modify the desired fields.

Field	Instruction
Sample Name	Edit the unique identifier for the sample.
Description	Enter or edit a preferred description of the sample.
Sample Type	Select a different sample type:  • Sample—A forensic or reference sample  • ForenSeq Positive Amplification Control—The control DNA (NA24385, or HL60)  • Negative Amplification Control—The negative control (water)  • Reagent Blank—The extraction control
Assay Type*	Select a different library prep method:  • ForenSeq mtDNA Whole Genome—Libraries prepared with the ForenSeq mtDNA Whole Genome Kit  • ForenSeq mtDNA Control Region—Libraries prepared with the ForenSeq mtDNA Control Region Kit  • mtDNA Custom—Libraries prepared with a custom or third-party assay  • ForenSeq Kintelligence—Libraries prepared with a ForenSeq Kintelligence Kit  • ForenSeq MainstAY—Libraries prepared with a ForenSeq MainstAY Kit
Index 1*	Select the name of a different i7 index adapter.
Index 2*	Select the name of a different i5 index adapter.
Index Combination*	Select the name of a different premixed i7 and i5 index adapter.
Run Analysis Method*	Select a different run analysis method.
Initial Projects*	Enter additional projects or remove current projects. For detailed instructions, see <i>Modify Project Assignments</i> .
Sample Plate	Edit or enter a preferred name for the 96-well plate.
Sample Well Position	Edit or enter the position of the well the sample occupies.

<sup>\*</sup> Editing these fields requires a run state of Created. To edit them after run completion, remove the sample from the run, make the necessary edits, and add the sample to the run again.

4. Select **Save** to apply the changes.

#### Modify Project Assignments

When the run state is Created, you can change which projects a sample is assigned to or assign the sample to additional projects.

- 1. On the navigation rail, select **Runs**, and then select a run.
- 2. In the Actions column, select Edit Sample Details for the applicable sample.
- In the Edit Sample dialog box, edit the Initial Projects field as follows.
  - To remove a sample from a project, select X next to the project name.
  - To add a sample to a project, enter the project name.
  - To add a sample to multiple projects, enter multiple project names.
  - To create a project for the sample, enter a new project name.
- 4. Select Save.

#### Remove Samples from a Run

- 1. On the navigation rail, select **Runs**, and then select a run.
- 2. Select the checkbox of each sample you want to remove from the run.
- 3. Select Remove Samples.
- 4. In the Remove Samples dialog box, select Confirm.

#### Monitor a Run

During a run, the Run Details page indicates which cycle the run is on and displays metrics so you can monitor run status and performance. Green indicates passing metrics and orange indicates metrics that require further evaluation. These metrics are a subset of the run metrics used to evaluate performance after sequencing is complete.

- 1. On the navigation rail, select **Runs**, and then select a run.
- 2. Select Quality Metrics.
  - Review the values for cluster density and clusters passing filter.
  - Review the phasing and prephasing values for Read 1.
  - Check Read 1, Index 1, Index 2, and Read 2 to confirm that all four reads passed.
  - Check Overall Intensity and 0 Discordant Loci to confirm that the HSC passed.

For more information, see Quality Metrics on page 68.

- 3. Select Close to return to the run.
- 4. Select Sample Representation.
  - Compare the samples to the sample read count guideline to confirm performance.
  - Review the number of reads in parentheses next to the sample name or hover over each sample.

For more information, see *Sample Representation* on page 75.

5. Select **Close** to return to the run.

# **Data Analysis**

Review Run Metrics	30
Review mtDNA Results	31
Review Kintelligence Results	33
Review MainstAY Results	
Modify Kintelligence and MainstAY Results	
Add Analyzed Samples to a Project	37
Reanalyze Samples	37

#### **Review Run Metrics**

Run metrics provide run quality, results for positive and negative controls, and sample representation for an overall view of run performance. Green indicates passing metrics and orange indicates metrics that require further evaluation.

For help with metrics requiring further evaluation, see *Troubleshoot a Run* on page 63. For detailed descriptions of each metric, see *Run Metrics* on page 68.

- 1. On the navigation rail, select **Projects**, and then select a project.
- 2. Select QPN for a run in the project.
- 3. On the Quality Metrics page, do as follows.
  - Review the values for cluster density and clusters passing filter.
  - Review the phasing and prephasing values for Read 1.
  - Check Read 1, Index 1, Index 2, and Read 2 to confirm that all four reads passed.
  - Check Overall Intensity and 0 Discordant Loci to confirm that the HSC passed.
- 4. Select Positive & Negative Controls.
  - Under ForenSeq Positive Amplification Controls, review the number of positions called or SNPs or STRs typed, discordant loci, and pass filter reads for each positive control.
    - Expand **Discordance** to see the Discordance table.
    - Use the scroll arrows to move through all positions in the table.
  - Under Reagent Blanks, review the number of positions called, SNPs typed, or STRs typed for each reagent blank.
  - Under Negative Amplification Controls, review the number of positions called, SNPs typed, or STRs typed for each negative control.
  - Select a sample name to go to the Results page.
- Select Sample Representation.
  - Compare the samples to the sample read count guideline to confirm performance.
  - Review the number of reads in parentheses next to the sample name or hover over each sample.
- 6. If any metrics do not pass, contact Verogen Technical Support. Collect the following items from D:\\llumina\\MiSeqFGxOutput to help Verogen with troubleshooting:

- InterOp folder
- runInfo.xml file
- runParameters.xml file

#### Review mtDNA Results

The Results page for each mtDNA sample includes filters, fields, and buttons for a dynamic view of results for the sample. The following sections describe how to use these features to review the results. For feature descriptions, see *Results Page for mtDNA* on page 13.

#### Filter Results

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Select Filter, and then select at least one of the following checkboxes:
  - Control Region—Zooms in on the control region of the mtGenome.
  - Interpretation Threshold-Indicates positions with reads between the AT and IT.
  - Variants—Indicates positions with base calls that differ from the reference.
  - Call—For reagent blanks and negative amplification controls, indicates positions with no calls to identify areas without coverage.
  - No call—For samples and positive amplification controls, indicates positions with calls to identify possible contamination.

The sample type determines whether Call or No Call filter is visible.

Reselect Filter to apply the selected filters.

The mtDNA navigator, position viewer, and coverage plot refresh to reflect the filtered data.

- 5. To clear the filters:
  - a. Select Active Filter.
  - b. Deselect the checkboxes.
  - c. Reselect Active Filter.

# Set Maximum Coverage

Setting a maximum coverage filters sample data by total read count.

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
- 3. In the Max Coverage field, enter a read number without commas.

The mtDNA navigator and coverage plot adjust to the specified value.

# Compare Samples

You can compare and filter results among 2-9 mtDNA samples in a project. During a comparison, the mtDNA navigator and coverage plot show data for the sample of interest in gray and a selected comparison sample in blue. The position viewer displays sequences and data for the sample of interest and selected comparison sample. The name of the sample of interest appears blue.

- 1. Navigate to the Results page:
  - a. On the navigation rail, select Projects.
  - b. Select a project, and then select a sample to be the sample of interest.
- 2. On the Samples page, select a sample.
  - This sample is the sample of interest.
- 3. On the Results page, select Compare.
- 4. In the Sample Compare dialog box, select the checkboxes of up to eight samples.
- 5. Select Submit.

The position viewer lists the selected comparison samples. The mtDNA navigator and coverage plot refresh to show data for the sample of interest and the first selected comparison sample. A green checkmark next to a selected comparison sample indicates no differences with the sample of interest.

- 6. If you are comparing more than two samples, select a comparison sample in the position viewer.
  - The mtDNA navigator and coverage plot refresh to compare the selected comparison sample to the sample of interest.
- 7. [Optional] Select Filter, and then select at least one of the following checkboxes:
  - Control Region—Zooms in on the control region of the mtGenome.
  - Interpretation Threshold—Indicates positions with reads between the AT and IT.
  - Variants—Indicates positions with base calls that differ from the reference.
  - Call—For reagent blanks and negative amplification controls, indicates positions with no calls to identify areas without coverage.
  - No call—For samples and positive amplification controls, indicates positions with calls to identify possible contamination.
  - Differences—Indicates call differences between the sample of interest and the selected comparison samples.

The sample of interest determines whether the Call or No Call filter is displayed. For example, when the sample of interest is a negative amplification control or reagent blank, the software displays the Call filter. When the sample of interest is a sample or positive amplification control, the software displays the No Call filter.

8. When you are done comparing samples, select **Stop Compare**.

#### **Edit Base Calls**

1. On the navigation rail, select **Projects**.

- 2. Select a project, and then select an mtDNA sample.
- 3. In the position viewer, select a position by using the scroll arrows or entering the position.
- Select the Call arrow, and then select a call. Scroll to see all available calls.
   The position viewer refreshes to reflect the modified call. An asterisk (\*) marks the modified call and the original call.

The position viewer refreshes to reflect the modified call. An asterisk (\*) marks the modified call and the original call appears below the Call list.

# Review Kintelligence Results

The Results page for each Kintelligence sample includes filters and sortable categories to organize and locate loci for review. The following sections describe how to use these features to review results.

- For information on the Results page, see Results Page for Kintelligence on page 16.
- For information on how the software estimates gender and contributor status, see Supporting Information on page 85.

By default, the Results page displays results with the Typed filter and all SNP Type filters applied except X-SNP and Y-SNP. Sort 1 is set to descending Allele Count, Sort 2 is set to ascending Intralocus Balance, and Sort 3 is set to ascending Intensity.

#### Filter Loci

- 1. On the navigation rail, select Projects.
- 2. Select a project, and then select a sample.

The Results page appears.

3. Under Sample Details, select any of the following filters to view only the desired loci:

Filter	Effect
Chromosome	Display loci located on the selected chromosome only.
Typed	Display only loci that are typed.
Untyped	Display only loci that are untyped.
Homozygous	Display only loci that have the same allele.
Heterozygous	Display only loci that have different alleles.
QC Indicators	Displays only loci flagged with the selected QC indicators.
SNP Type	Displays only loci that belong to the selected SNP types.

- To remove filters, select the following options:
  - Clear Selection—Clear all filters in a section.
  - Default View-Clear selected filters and restore the default filters.

#### Sort Loci

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.

The Results page appears.

3. Under Sample Details, select the Sort 1, Sort 2, or Sort 3 list, and then select one of the following options:

Effect
Do not apply any sorting.
Sort loci by intensity value.
Sort loci by intralocus balance percentage.
Sort loci by amplicon size in bp.
Sort loci by chromosome number or letter.
Sort loci by SNP category.
Sort loci by rsID.
Sort loci by the number of alleles.

- To further sort loci into ascending or descending order, select the arrow next to each sort list.
- 5. To clear your sorting selections and restore the defaults, select **Default View**.

#### Find a Locus

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select a sample.
  - The Results page appears.
- In the Search field, enter the whole or partial name of the locus you want to find, and then press Enter.The Results page displays loci matching the search term.

#### Review MainstAY Results

The Results page for each MainstAY sample includes filters and sorting to organize and locate loci for review. By default, the histogram view sorts loci in ascending order by chromosome without any filters applied. The following sections describe how to use these features to review the results. For feature descriptions, see *Results Page for MainstAY* on page 17.

#### Filter Loci

- 1. On the navigation rail, select Projects.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- In the Histogram View, on the filter pane, select All, Autosomal STRs, or Y STRs.
- 4. Under the QC Indicators section, select the check box associated with the filter that needs to be applied.
  - To apply all filters except No QC Indicators, select Select All.
  - To clear all filters, select Clear Selection.

The software displays only the loci flagged with the selected indicators. For indicator descriptions, see *Locus Quality Indicators* on page 18.

#### Sort Loci

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Under Histogram View, select All, Autosomal STRs, or Y STRs.
- 4. Select **Sort By**, and then select one of the following options:

Sort Option	Effect
Allele Count	Sort loci by the number of alleles.
Amplicon Size	Sort loci by amplicon size in bp.
Chromosome	Sort loci by chromosome number or letter.
Intralocus Balance	Sort loci by intralocus balance percentage.
Intensity	Sort loci by intensity value.
STR Name	Sort loci by locus name.
STR Type	Sort loci by aSTRs and Y-STRs.

5. To further sort loci into ascending or descending order, select the arrow to the left of the Sort By list.

#### Find a Locus

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Under Histogram View, select All, Autosomal STRs, or Y STRs.
- 4. In the Search field, enter the whole or partial name of the locus you want to find.
  - The Results page displays loci matching the search term.
- 5. To return to the full view, clear the Search field.

# Modify Kintelligence and MainstAY Results

The software allows you to modify Kintelligence and MainstAY results by typing and untyping alleles, changing the sample gender, and changing the contributor status. For MainstAY results, you can add comments to loci explaining any changes. For instructions, see *Comment on a MainstAY Locus* on page 54.

A green toggle key and green histogram bar indicate a typed allele. A gray toggle key and bar indicate an untyped allele.

# Type or Untype a Kintelligence Allele

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a Kintelligence sample.
- 3. On the Results page, under Sample Details, select a locus.

The locus details appear.

- 4. Type or untype an allele using one of the following methods:
  - · Under Typed, select a toggle key.
  - · On the histogram, select a bar.

The User Modified icon appears, indicating the modification. Reverting the call removes the icon.

#### Type or Untype an MainstAY Allele

- 1. On the navigation rail, select Projects.
- 2. Select a project, and then select an STR sample to open the Results page.
- 3. Under Histogram View, select All, Autosomal STRs, or Y STRs.
- 4. Select the locus you want to edit.

The locus details appear.

- 5. Type or untype an allele using one of the following methods:
  - Under Typed, select a toggle key.
  - On the histogram, select a bar.

The User Modified icon appears, indicating the modification. Reverting the call removes the icon.

# Change the Contributor Status

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a Kintelligence or MainstAY sample.

The Results page appears.

- 3. In the Sample Overview section, select the Contributor Status list, and then select a status:
  - Single-Source—The sample is from one contributor.
  - Mixture—The sample is from multiple contributors.
  - Inconclusive—Whether the sample is from one contributor or multiple contributors is unknown.

# Change the Sample Gender

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a Kintelligence or MainstAY sample.

The Results page appears

- 3. In the Sample Overview section, select the Gender list, and then select the applicable chromosomes:
  - XX–The sample is from a female contributor.

- XY—The sample is from a male contributor.
- Inconclusive—The biological sex of the contributor is unknown.

# Add Analyzed Samples to a Project

You can add analyzed samples in a project to other projects. After adding a sample to another project, sample modifications are *not* applied in the newly assigned project.

- 1. On the navigation rail, select **Projects**, and then select a project.
- 2. On the Samples page, select **Add Samples**.
- 3. In the Add Samples dialog box, in the **Search** field, enter a sample name.
- 4. In the search results, select the checkbox of each sample you want to add to the project. Each selection appears in a summary list with the corresponding analysis method.
- 5. To remove a sample, deselect the checkbox or select **X** in the summary.
- 6. Select Add to add the samples.

# Reanalyze Samples

Reanalyze samples to apply an analysis method with user-defined parameters. The software preserves results from all analysis methods.

- 1. On the navigation rail, select Projects, and then select a project.
- 2. On the Samples page, select the checkbox of each sample you want to reanalyze.
- 3. Select Reanalyze.
- 4. In the Reanalyze Samples dialog box, choose one of the following options:

Reanalysis Option	Instruction
Assign the same analysis method to all samples selected for reanalysis.	Select a method from the Analysis Methods list, and then select <b>Apply to All</b> .
Assign an analysis method to each sample selected for reanalysis.	Select a method from each list in the New Analysis Method column.

The software only lists analysis methods that are compatible with the assay used to prepare the sample.

- Select Next.
- 6. Review the sample names, previous analysis methods, and new analysis methods.
- 7. To make changes, select **Previous**.
- 8. Select **Save** to initiate reanalysis.

In the Projects sidebar, an icon appears next to each sample undergoing reanalysis.

9. When reanalysis is complete, select a sample to view results.

The Results page displays results for the most recent analysis.

## Switch Analysis Methods

A Results page shows results for the analysis method displayed in the Select an Analysis Method list. For reanalyzed samples, the list includes all analysis methods applied to the sample with the most recent first. Select different analysis methods to see other results.

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. In the Select an Analysis Method list, select an analysis method.
  - The Results page is refreshed to show results for the selected method.

# Report Generation

Generating mtDNA Reports	. 39
Generating Kintelligence Reports	. 43
Generating MainstAY Reports	45

# Generating mtDNA Reports

The software compiles mtDNA results into the following reports:

- CODIS—Variants from the selected samples for upload to CODIS.
- **EMPOP**—Specimen data, originating agency identifiers (ORIs), and other DNA profile information for upoad to EMPOP.
- Project—Results for all samples in a project.
- Sample—Results for one sample.

The following sections describe how to generate each type of report. For report contents, see *File Formats for mtDNA Reports* on page 78.

## Create a CODIS Report

Create a CODIS report for one sample or multiple samples. Creating a report for one sample lets you specify extra specimen information and choose which positions to include in the report. When creating a report for multiple samples, the software automatically includes all positions for each sample.

#### Create a CODIS Report for One Sample

- 1. On the navigation rail, select Projects.
- 2. Select a project, and then select a sample.
- 3. Select **Reports**, and then select **mtDNA CODIS Report**.
  - The Create CODIS Report dialog box appears.
- [Optional] Select CODIS Report Defaults to edit default report settings.
  - For more information, see *CODIS Report Defaults* on page 60.
- 5. Edit the following CMF header fields or accept the default values.
  - Destination Laboratory ORI-Identifier for the agency that processes the CODIS report.
  - Source Laboratory ORI-Identifier for the agency that processes the specimen.
  - Submission User ID—Identification of the person who generates the CODIS report.
  - Batch ID—Identifier that tracks where a DNA profile originated.
- 6. To include the name of the library prep kit in the report, select Include Kit Name in Report.
- 7. Select Next.
- 8. In the Specimen ID field, enter a unique identifier for the specimen or accept the default sample name.
- 9. In the Specimen Category list, select a preconfigured category or select **Custom** and enter the category name.

#### 10. [Optional] Select More Options.

- a. Complete the following fields:
  - Source ID—Specify whether the identity of the specimen contributor is known.
  - NCIC Number-Enter the unique number for storage in the NCIC system.
  - ViCAP number—Enter the unique number for storage in the ViCAP system.
  - Case ID—Enter the unique identifier for the case.
  - Case ID—Enter the law enforcement identifier associated with the specimen.
  - Specimen Comment—Enter any additional information about the specimen.
- b. If the specimen has possible allele dropout or other conditions requiring additional information to help with match resolution, select **Partial Profile**.
- Select Next.
- 12. Select the toggle key for each position you want to remove from the report.
  - By default, all positions are included.
- 13. Select Create to generate the report.
- 14. Select **Reports** to open the Reports page.

The Reports page lists all reports for the project with the new report on top.

15. In the Actions column, select Download Report to download the report as a .cmf file.

#### Create a CODIS Report for Multiple Samples

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select Reports.
- Select Create Report, and then select mtDNA CODIS Report.

The Create CODIS Report dialog box appears.

- 4. Select the checkbox of each sample you want to include in the report.
  - Scroll and use the arrows to move through the complete list of samples in the project.
  - Under Summary, select Delete to remove samples from the report.
- Select Next.
- [Optional] Select CODIS Report Defaults to edit default report settings.

For more information, see *CODIS Report Defaults* on page 60.

- 7. Edit the following CMF header fields or accept the default values:
  - Destination Laboratory ORI—Identifier for the agency that processes the CODIS report.
  - Source Laboratory ORI—Identifier for the agency that processes the specimen.
  - Submission User ID—Identification of the person who generates the CODIS report.
  - Batch ID—Identifier that tracks where a DNA profile originated.
- 8. To include the name of the library prep kit in the report, select Include Kit Name in Report.
- Select Next.

- 10. In the Specimen ID field for each sample, enter a unique identifier for the specimen cited in the report.
- 11. Assign specimen categories per the following table. If you select Custom, enter a name for the custom specimen category.

Option	Instruction
Assign the same category to all samples.	In the Specimen Category list at the top of the dialog box, select a category, and then select <b>Apply to All Specimens</b> .
Assign a different category to each sample.	In the Specimen Category list for each sample, select the applicable category.

12. Select Create to generate the report.

The software returns to the Reports page, where the new report appears at the top of the list.

13. In the Actions column, select **Download Report** to download the report as a .cmf file.

#### Create an EMPOP Report

Create an EMPOP report for one or more samples in a project. When creating a report for one sample, you can start from the Results or Reports page.

#### Create an EMPOP Report for One Sample from Results

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- Select Reports, and then select mtDNA EMPOP Report to generate the report.
- 4. Select **Reports** to open the Reports page.
  - The Reports page lists the new report on top.
- 5. In the Actions column, select **Download** to download the report as a .txt file.

#### Create an EMPOP Report for One Sample from Reports

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select Reports.
  - The Reports page appears.
- 3. Select Create Report, and then select mtDNA EMPOP Report.
- 4. In the Create mtDNA EMPOP Report dialog box, select the checkbox of the applicable sample.
- 5. Select **Confirm** to generate the report.
  - The returns to the Reports page, where the new report appears at the top of the list.
- 6. In the Actions column, select **Download Report** to download the report as a .txt file.

## Create an EMPOP Report for Multiple Samples

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select Reports.

The Reports page appears.

- Select Create Report, and then select mtDNA EMPOP Report.
- 4. In the Create EMPOP Report dialog box, select the checkbox of each sample to include in the report.
  - Scroll and use the arrows to move through the complete list of samples in the project.
  - Under Summary, select Delete to remove samples from the report.
- 5. Select **Confirm** to generate the report.

The returns to the Reports page, where the new report appears at the top of the list.

6. In the Actions column, select **Download Report** to download the report as a .txt file.

#### Create a Project Report

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select Reports.
- On the Reports page, select Create Report, and then select Project.
   If the software is analyzing samples, project report creation is disabled and a notification is displayed.
- 4. When prompted, select **Yes** to create the report.

The software generates the report and returns to the Reports page, where the new report appears at the top of the list.

If a project includes samples with different assay types, the software generates one project report per assay.

5. In the Actions column, select **Download** to download a report as an .xmlx file.

## Create Sample Reports

Create one sample report at a time or simultaneously generate multiple sample reports, one for each selected sample in a project. When creating a report for one sample, you can start from the Results or Reports page.

#### Create One Sample Report from Results

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Select Reports, and then select Sample Report to generate the report.
- 4. Select Reports to open the Reports page.
  - The Reports page lists the new report on top.
- 5. In the Actions column, select **Download** to download the report as an .xmlx file.

#### Create One Sample Report from Reports

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select Reports.

The Reports page appears.

- 3. Select Create Report, and then select Sample.
- 4. In the Create Sample Report dialog box, select the checkbox of the sample you want to create a report for. Scroll and use the arrows to see all samples.
- 5. Select Confirm to generate the report.

The software returns to the Reports page, where the new report appears at the top of the list.

6. In the Actions column, select **Download** to download a report as an .xmlx file.

#### Create Multiple Sample Reports

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select Reports.

The Reports page appears.

- Select Create Report, and then select Sample.
- 4. In the Create Sample Report dialog box, select the checkbox of each sample you want to create a report for.
  - Scroll and use the arrows to move through the complete list of samples in the project.
  - · Use the Delete icons under Summary to remove samples from the report.
- Select Confirm to generate the reports.

The software returns to the Reports page, where the new reports appear at the top of the list.

6. In the Actions column, select **Download** to download a report as an .xmlx file.

# Generating Kintelligence Reports

The software compiles Kintelligence results and other information into the following reports:

- Phenotype and Ancestry

  —Estimates of hair color, eye color, and biogeographical ancestry based on aSNP and pSNP data.
- GEDmatch PRO—SNP genotype calls reported on the plus strand for upload to GEDmatch PRO.
- Project—Results for all samples in a project.
- Sample—Results for one sample.

The following sections describe how to generate each type of report. For report contents, see *File Formats for Kintelligence Reports* on page 80.

## Create a Phenotype and Ancestry Report

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select a sample.

The Results page appears.

- 3. Select Create Report, and then select Phenotype & Ancestry to generate the report.
- Select Reports to open the Reports page.

The Reports page lists the new report on top.

5. In the Actions column, select **Download** to download the report as an .xmlx file.

## Create a GEDmatch PRO Report

- 1. On the navigation rail, select Projects.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Select Reports, and then select GEDmatch PRO to generate the report.
- 4. Select **Reports** to open the Reports page.
  - The Reports page lists the new report on top.
- 5. In the Actions column, select **Download** to download the report as a .txt file.

## Create a Project Report

- 1. On the navigation rail, select Projects.
- 2. Select a project, and then select Reports.
- 3. On the Reports page, select Create Report, and then select Project.
  - If the software is analyzing samples, project report creation is disabled and a notification is displayed.
- 4. When prompted, select **Yes** to create the report.
  - The software generates the report and returns to the Reports page, where the new report appears at the top of the list

If a project includes samples with different assay types, the software generates one project report per assay.

5. In the Actions column, select **Download** to download a report as an .xmlx file.

## Create Sample Reports

Create one sample report at a time or simultaneously generate multiple sample reports, one for each selected sample in a project.

#### Create One Sample Report

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Select Reports, and then select Sample Report to generate the report.
- 4. Select Reports to open the Reports page.
  - The Reports page lists the new report on top.
- 5. In the Actions column, select **Download** to download the report as an .xmlx file.

#### Create Multiple Sample Reports

1. On the navigation rail, select **Projects**.

Select a project, and then select Reports.

The Reports page appears.

- 3. Select Create Report, and then select Sample.
- 4. In the Create Sample Report dialog box, select the checkbox of each sample you want to create a report for.
  - Scroll and use the arrows to move through the complete list of samples in the project.
  - Use the Delete icons under Summary to remove samples from the report.
- 5. Select **Confirm** to generate the reports.

The software returns to the Reports page, where the new reports appear at the top of the list.

6. In the Actions column, select **Download** to download a report as an .xmlx file.

# Generating MainstAY Reports

The software compiles STR results into the following reports:

- CODIS-Variants from the selected samples for upload to CODIS.
- Project—Results for all samples in a project.
- Sample—Results for one sample.

The following sections describe how to generate each type of report. For report contents, see *File Formats for MainstAY Reports* on page 83.

#### Create a CODIS Report

Create a CODIS report for one sample or multiple samples. Creating a report for one sample lets you specify extra specimen information and choose which positions to include in the report. When creating a report for multiple samples, the software automatically includes all positions for each sample.

#### Create a CODIS Report for One Sample

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select a sample.
- 3. Select Reports, and then select CODIS Report.

The Create CODIS Report dialog box appears.

4. [Optional] Select CODIS Report Defaults to edit default report settings.

For more information, see CODIS Report Defaults on page 60.

- 5. Edit the following CMF header fields or accept the default values.
  - Destination Laboratory ORI

    —Identifier for the agency that processes the CODIS report.
  - Source Laboratory ORI-Identifier for the agency that processes the specimen.
  - Submission User ID—Identification of the person who generates the CODIS report.
  - Batch ID—Identifier that tracks where a DNA profile originated.
- To include the name of the library prep kit in the report, select Include Kit Name in Report.
- Select Next.

- 8. In the Specimen ID field, enter a unique identifier for the specimen or accept the default sample name.
- 9. In the Specimen Category list, select a preconfigured category or select Custom and enter the category name.
- 10. [Optional] Select More Options.
  - a. Complete the following fields:
    - Source ID—Specify whether the identity of the specimen contributor is known.
    - NCIC Number—Enter the unique number for storage in the NCIC system.
    - ViCAP number—Enter the unique number for storage in the ViCAP system.
    - Case ID—Enter the unique identifier for the case.
    - Case ID—Enter the law enforcement identifier associated with the specimen.
    - Specimen Comment—Enter any additional information about the specimen.
  - b. If the specimen has possible allele dropout or other conditions requiring additional information to help with match resolution, select **Partial Profile**.
- 11. Select Next.
- 12. Select the toggle key for each position you want to remove from the report.
  - By default, all positions are included.
- 13. Select Create to generate the report.
- 14. Select Reports to open the Reports page.

The Reports page lists all reports for the project with the new report on top.

15. In the Actions column, select Download Report to download the report as a .cmf file.

#### Create a CODIS Report for Multiple Samples

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select Reports.
- 3. Select Create Report, and then select CODIS Report.

The Create CODIS Report dialog box appears.

- 4. Select the checkbox of each sample you want to include in the report.
  - Scroll and use the arrows to move through the complete list of samples in the project.
  - Under Summary, select Delete to remove samples from the report.
- Select Next.
- 6. [Optional] Select CODIS Report Defaults to edit default report settings.

For more information, see *CODIS Report Defaults* on page 60.

- 7. Edit the following CMF header fields or accept the default values:
  - Destination Laboratory ORI—Identifier for the agency that processes the CODIS report.
  - Source Laboratory ORI-Identifier for the agency that processes the specimen.
  - Submission User ID—Identification of the person who generates the CODIS report.
  - Batch ID—Identifier that tracks where a DNA profile originated.

- 8. To include the name of the library prep kit in the report, select Include Kit Name in Report.
- Select Next.
- 10. In the Specimen ID field for each sample, enter a unique identifier for the specimen cited in the report.
- 11. Assign specimen categories per the following table. If you select Custom, enter a name for the custom specimen category.

Option	Instruction
Assign the same category to all samples.	In the Specimen Category list at the top of the dialog box, select a category, and then select <b>Apply to All Specimens</b> .
Assign a different category to each sample.	In the Specimen Category list for each sample, select the applicable category.

12. Select Create to generate the report.

The software returns to the Reports page, where the new report appears at the top of the list.

13. In the Actions column, select **Download Report** to download the report as a .cmf file.

#### Create a Project Report

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select Reports.
- On the Reports page, select Create Report, and then select Project.
   If the software is analyzing samples, project report creation is disabled and a notification is displayed.
- 4. When prompted, select **Yes** to create the report.

The software generates the report and returns to the Reports page, where the new report appears at the top of the list.

If a project includes samples with different assay types, the software generates one project report per assay.

5. In the Actions column, select **Download** to download a report as an .xmlx file.

## Create Sample Reports

Create one sample report at a time or simultaneously generate multiple sample reports, one for each selected sample in a project.

#### Create One Sample Report

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Select Reports, and then select Sample Report to generate the report.
- 4. Select Reports to open the Reports page.

The Reports page lists the new report on top.

5. In the Actions column, select **Download** to download the report as an .xmlx file.

## Create Multiple Sample Reports

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select Reports.
  - The Reports page appears.
- 3. Select Create Report, and then select Sample.
- 4. In the Create Sample Report dialog box, select the checkbox of each sample you want to create a report for.
  - · Scroll and use the arrows to move through the complete list of samples in the project.
  - Use the Delete icons under Summary to remove samples from the report.
- 5. Select **Confirm** to generate the reports.

The software returns to the Reports page, where the new reports appear at the top of the list.

6. In the Actions column, select **Download** to download a report as an .xmlx file.

# Run, Project, and Sample Management

Create a Project	49
Edit a Project	
Obtain the Analysis Folder Path	50
Printing	50
History	

# Create a Project

Create a project to group a selection of samples. For example: create a project of controls to compare controls from multiple projects, or group samples in a comparison to generate a project report of compared samples.

Although possible, combining samples with different assay types in a project is not recommended.

- 1. On the navigation rail, select **Projects**.
- 2. Select Create Project.
- 3. In the Create Project dialog box, complete the following fields:
  - Project Name—Enter a preferred name to identify the project.
  - Description—Optionally enter a description of the project.
- 4. Select Save to create the project.

The new project opens to the Samples page.

5. Select Add Samples.

You can add any samples in the software to the project, analyzed or not.

- 6. In the Add Samples dialog box, in the **Search** field, enter a sample name.
- 7. In the search results, select the checkbox of each sample you want to add to the project.

Each selection appears in a summary list with the corresponding analysis method.

- 8. To remove a sample, deselect the checkbox or select **X** in the summary.
- 9. Select Add to add the samples.

The Samples page lists the added samples. Runs the samples belong to appear in the Project sidebar.

# Edit a Project

- 1. On the navigation rail, select **Projects**, and then select a project.
- Select Edit.
- 3. In the Edit Project dialog box, do as follows.
  - Edit the project name.
  - Add or modify the optional project description.
- 4. Select Save to apply the changes.

# Obtain the Analysis Folder Path

The analysis folder resides on the ForenSeq Universal Analysis Server and contains the Variant Call Format (VCF) and Binary Alignment Map (BAM) files for the analysis.

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select a sample.
  - The Results page appears.
- 3. Select the Information icon.
- 4. The path is displayed under **Analysis Folder Path**.

## **Printing**

The Print icon converts the following pages into a printable format:

- Run Details—Print the header of a Run Details page, the table listing samples in the run, run quality metrics, and sample representation.
- · Samples—Print the Samples page of a project.
- Results—Applicable to Kintelligence samples only, print the header of a Results page and the Sample Overview section.

#### Print a Run Details Page

- 1. On the navigation rail, select **Runs**, and then select a run.
- On the Run Details page, select Print.
  - The software displays a printable view of the Run Details page.
- 3. Select **Print** to open the Print dialog box.
- 4. Choose printing options and print.
- Select Back to Run.

## Print a Samples Page

- 1. On the navigation rail, select **Projects**, and then select a project.
- On the Project sidebar, select Print.
  - The software displays a printable view of the Samples page.
- 3. Select **Print** to open the Print dialog box.
- 4. Choose printing options and print.
- 5. Select Back to Project.

#### Print a Sample Overview Section

- 1. On the navigation rail, select Projects.
- 2. Select a project, and then select a Kintelligence sample.

- On the Results page, select Print Sample Results.
  - The software displays a printable view of the Sample Overview section.
- 4. Select **Print** to open the Print dialog box.
- 5. Choose printing options and print.
- 6. Select Back to Project.

## History

When the history setting is enabled, you can review system events and user actions recoded for a run, project, or sample and add comments. For example, you can comment on a sample to explain a call modification. For instructions on enabling the history setting, see *System Settings* on page 61.

The software also records system events and user actions for loci in MainstAY samples and allows commenting. Loci history and commenting is always available, regardless of whether the history setting is enabled.

### Review Run Activity

- 1. On the navigation rail, select Runs.
- Select a run, and then select Activity at the top of the page.
- 3. In the Run Activity dialog box, review system events and user actions for the run.
  - The dialog box lists system events and user actions from oldest to newest.
- 4. Navigate using the scroll bar and arrows. Narrow results using the following fields:
  - Search—Search the current history.
  - Start Date and End Date—View the history within a specified time frame.
  - Show System Events—See system events and user actions for the history or see user actions only.
- 5. When you are finished, select Close.

## **Review Project Activity**

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select Project Activity at the top of the page.
- 3. In the Project Activity dialog box, review system events and user actions for the project.
  - The dialog box lists system events and user actions from oldest to newest.
- 4. Navigate using the scroll bar and arrows. Narrow results using the following fields:
  - Search—Search the current history.
  - Start Date and End Date—View the history within a specified time frame.
  - Show System Events—See system events and user actions for the history or see user actions only.
- 5. When you are finished, select Close.

## Review Sample Activity

1. On the navigation rail, select **Projects**.

- 2. Select a project, and then select a sample.
- 3. Select Sample Activity.
- 4. In the Sample Activity dialog box, review system events and user actions for the sample.
  - The dialog box lists system events and user actions from oldest to newest.
- 5. Navigate using the scroll bar and arrows. Narrow results using the following fields:
  - Search—Search the current history.
  - Start Date and End Date—View the history within a specified time frame.
  - Show System Events—See system events and user actions for the history or see user actions only.
- 6. When you are finished, select Close.

#### Review MainstAY Locus Activity

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select a MainstAY sample.
  - The Results page appears.
- 3. Under Histogram View, select All, Autosomal STRs, or Y STRs.
- 4. Select a card to open the locus details.
- 5. Select the User Actions tab.

The tab lists system events and user actions from oldest to newest.

- 6. Review system events and user actions for the locus.
- 7. Navigate using the arrows. Narrow results using the following fields:
  - Search—Search the current history.
  - Start Date and End Date—View the history within a specified time frame.
  - Show System Events—See system events and user actions or see user actions only.
- 8. When you are finished, select Close.

#### Comment on a Run

- 1. On the navigation rail, select **Runs**.
- 2. Select a run, and then select Activity at the top of the page.
  - The Run Activity dialog box appears.
- 3. To add a comment to the run history:
  - a. In the Add New Comment field, type a comment.
  - b. Select Add to save the comment.

The comment appears as an event at the top of the history.

- 4. To comment on an entry in the run history:
  - a. Select the system event or user action. Use the fields and scrolling to help locate the entry, if necessary.

- b. In the Edit Comment box, type a comment.
- c. Select Save to add the comment to the system event or user action.

The comment appears as the last item for the entry.

5. When you are finished, select Close.

## Comment on a Project

- 1. On the navigation rail, select **Projects**.
- 2. Select a project, and then select **Project Activity** at the top of the page.

The Project Activity dialog box appears.

- 3. To add a comment to the project history:
  - a. In the Add New Comment field, type a comment.
  - b. Select Add to save the comment.

The comment appears as an event at the top of the history.

- 4. To comment on an entry in the project history:
  - a. Select the system event or user action. Use the fields and scrolling to help locate the entry, if necessary.
  - b. In the Edit Comment box, type a comment.
  - c. Select **Save** to add the comment to the system event or user action.

The comment appears as the last item for the entry.

5. When you are finished, select Close.

## Comment on a Sample

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select a sample.
- 3. Select Sample Activity.

The Sample Activity dialog box appears.

- To add a comment to the sample history:
  - a. In the Add New Comment field, type a comment.
  - b. Select Add to save the comment.

The comment appears as an event at the top of the history.

- 5. To comment on an entry in the sample history:
  - a. Select the system event or user action. Use the fields and scrolling to help locate the entry, if necessary.
  - b. In the Edit Comment box, type a comment.
  - Select Save to add the comment to the system event or user action.

The comment appears as the last item for the entry.

6. When you are finished, select Close.

#### Comment on a MainstAY Locus

Add a comment to a locus to, for example, explain why you changed an allele from untyped to typed. You can also edit or delete existing comments.

- 1. On the navigation rail, select **Projects**.
- Select a project, and then select a MainstAY sample.The Results page appears.
- 3. Under Histogram View, select All, Autosomal STRs, or Y STRs.
- 4. Select a card to open the locus details.
- 5. Select the **User Actions** tab.
- 6. In the Add New Comment field, type a comment.
- 7. Select **Add** to add the comment to the locus.
  - To edit the comment, select it, modify it as needed in the Edit Comment field, and select Save.
  - To delete the comment, select it, and then select Delete Comment.

# Software Settings

Settings Overview	55
Analysis Methods	55
CODIS Report Defaults	60
User Management	
System Settings	61
Enable or Disable Dark Mode	

# **Settings Overview**

Use the following settings, which are available from the Settings icon on the navigation rail, to configure the software. All settings except Dark Mode are visible to administrators only.

- Analysis Methods—Create and manage analysis methods.
- CODIS Report Defaults—Set default values to prepopulate CODIS reports.
- User Management–Add users and assign permissions.
- System Settings—Set visibility into system events and user actions.
- Dark Mode—Change the color scheme of the user interface.

# **Analysis Methods**

An analysis method is a collection of settings that inform how the software analyzes sequencing data and performs variant and allele calling. The software ships with a default Verogen analysis method for each of the following assays.

Assay Type	Default Analysis Method
ForenSeq Kintelligence	Verogen Kintelligence Analysis Method
ForenSeq mtDNA Whole Genome	Verogen mtDNA Whole Genome Analysis Method
ForenSeq mtDNA Control Region	Verogen mtDNA Control Region Analysis Method
mtDNA Custom	Verogen mtDNA Custom Analysis Method
ForenSeq MainstAY	ForenSeq MainstAY Analysis Method

# Create an Analysis Method

The default Verogen analysis methods cannot be modified, but you can create analysis methods with the desired settings using one of two methods: copy an existing analysis method and modify the settings or start with a blank analysis method.

Performance is supported when using the default analysis methods. Perform internal validation studies as needed to develop settings for new analysis methods.

For information on how Verogen determined the default settings in the mtDNA analysis methods and help determining alternative settings, see *Analysis Settings in ForenSeq Universal Analysis Software v2.0 - Verogen analysis methods offer fast, flexible and accurate variant calling or mitochondrial DNA data (document # VD202004).* For Kintelligence, see *SNP Typing in Universal Analysis Software and Kinship Estimation with GEDmatch PRO Technical Note (document # VD2020058) For MainstAY, see Robust thresholds and filters for reliable implementation of the ForenSeq MainstAY kit. (document # VD2021030).* 

#### Copy an Analysis Method

- 1. On the navigation rail, select **Settings**.
- 2. Select Analysis Methods.
- 3. Select Create Analysis Method.

The New Analysis Method dialog box appears.

- 4. In the New Analysis Method Options list, select Copy an Existing Analysis Method.
- 5. In the Source Analysis Method list, select an analysis method to use as a template.
- 6. In the New Analysis Method Name field, enter a unique name for the new analysis method.
- Select Save to create the analysis method.

The new method appears with the default settings from the copied method.

8. Modify the default settings as desired.

For instructions, see *Modify Analysis Settings* on page 56.

#### Use a Blank Analysis Method

- 1. On the navigation rail, select **Settings**.
- Select Analysis Methods.
- Select Create Analysis Method.

The New Analysis Method dialog box appears.

- 4. In the New Analysis Method Options list, select New Blank Analysis Method.
- 5. In the New Analysis Method Name field, enter a unique name for the new analysis method.
- 6. Select **Save** to create the analysis method.

The new method appears with blank settings.

7. Enter the desired settings for the analysis method.

For instructions, see *Modify Analysis Settings*.

#### Modify Analysis Settings

You can modify the following analysis settings, depending on the analysis method:

 For a user-created mtDNA analysis method, you can edit the AT, IT, minimum Q-score, and minimum read count. If an analysis method is copied from the Verogen mtDNA Custom Analysis Method, you can also modify the library type and amplicons.

- For a user-created Kintelligence analysis method, you can edit the AT, IT, and intralocus balance and specify which loci to include in the analysis.
- For a user-created MainstAY analysis method, you can edit the AT, IT, stutter filter, and intralocus balance and specify which loci to include in the analysis.

#### Edit an mtDNA Whole Genome or Control Region Analysis Method

- 1. If necessary, navigate to the analysis method you want to edit:
  - a. On the navigation rail, select **Settings**.
  - b. Select Analysis Methods.
  - c. In the Select an Assay Type list, select the assay that corresponds to the applicable analysis method.
  - d. Select the applicable analysis method.
- 2. To rename the analysis method, do as follows.
  - a. Select **Edit** to open the Edit Analysis Info dialog box.
  - b. In the Analysis Method Name field, enter a new name.
  - c. Select **Save** to apply the name change and close the dialog box.
- 3. Edit the following settings as desired:
  - In the Analytical Threshold and Interpretation Threshold boxes, adjust percentages.
  - · In the Minimum Quality Score and Minimum Read Count boxes, adjust values.

All percentages and values are greater than or equal to the value entered.

Select Save to apply the new settings.

#### Edit an mtDNA Custom Analysis Method

- 1. If necessary, navigate to the analysis method you want to edit:
  - a. On the navigation rail, select Settings.
  - b. Select Analysis Methods.
  - c. In the Select an Assay Type list, select the assay that corresponds to the applicable analysis method.
  - d. Select the applicable analysis method.
- 2. To rename the analysis method, do as follows.
  - a. Select Edit to open the Edit Analysis Info dialog box.
  - b. In the Analysis Method Name field, enter a new name.
  - c. Select **Save** to apply the name change and close the dialog box.
- Edit the following settings as desired:
  - In the Analytical Threshold and Interpretation Threshold boxes, adjust percentages.
  - In the Minimum Quality Score and Minimum Read Count boxes, adjust values.

All percentages and values are greater than or equal to the value entered.

- 4. In the Library Type list, select one of the following options:
  - Non-Directed Sequencing—Either strand can be sequenced as Read 1 or Read 2. This option is the default.
  - Primer-Directed Sequencing—One strand is sequenced as Read 1 and the opposite strand as Read 2.

The software lists amplicons available for the selected library type.

- 5. Modify the amplicons:
  - To remove an amplicon, in the Actions column, select Delete.
  - To edit an amplicon, in the Actions column, select Edit. In the Edit Amplicon dialog box, edit the desired fields, and then select Save.
  - To add an amplicon, select Add Amplicon. In the Add Amplicon dialog box, complete the necessary fields, and then select Confirm.

The user interface displays required fields with an asterisk (\*).

6. Select **Save** to apply the new settings.

#### Edit a Kintelligence Analysis Method

- 1. If necessary, navigate to the analysis method you want to edit:
  - a. On the navigation rail, select **Settings**.
  - b. Select Analysis Methods.
  - c. In the Select an Assay Type list, select the assay that corresponds to the applicable analysis method.
  - d. Select the applicable analysis method.
- To rename the analysis method, do as follows.
  - a. Select **Edit** to open the Edit Analysis Info dialog box.
  - b. In the Analysis Method Name field, enter a new name.
  - c. Select **Save** to apply the name change and close the dialog box.
- 3. Select the checkbox of at least one locus to edit.
  - Enter a whole or partial locus name in the Search Locus Name field to find a locus.
  - Select the checkbox of applicable SNP types to filter the list of loci.
  - · Select any column heading to sort the column.
  - Scroll to see all loci.

The software applies any new settings to all selected loci. If you want to apply different settings to different loci, select on locus at a time.

Select Edit Selected Loci.

The Edit Analysis Method dialog box appears.

5. To add or remove the loci from the analysis, select **Include**.

Green indicates that the loci are included. Gray indicates exclusion.

6. To modify the analysis thresholds, enter a value **0-100** in the applicable fields: Analytical Threshold, Interpretation Threshold, and Intralocus Balance.

The IT value must be greater than or equal to the AT value.

7. Select **Save** to apply the new settings.

#### Edit a MainstAY Analysis Method

- 1. If necessary, navigate to the analysis method you want to edit:
  - a. On the navigation rail, select **Settings**.
  - b. Select **Analysis Methods**.
  - c. In the Select an Assay Type list, select the assay that corresponds to the applicable analysis method.
  - d. Select the applicable analysis method.
- 2. To rename the analysis method, do as follows.
  - a. Select Edit to open the Edit Analysis Info dialog box.
  - b. In the Analysis Method Name field, enter a new name.
  - c. Select **Save** to apply the name change and close the dialog box.
- 3. Select the checkbox of at least one locus to edit.
  - Enter a whole or partial locus name in the Search Locus Name field to find a locus.
  - Select the checkbox of applicable STR types to filter the list of loci.
  - Select any column heading to sort the column.
  - Scroll to see all loci.

The software applies any new settings to all selected loci. If you want to apply different settings to different loci, select on locus at a time.

Select Edit Selected Loci.

The Edit Analysis Method dialog box appears.

- 5. To add or remove the loci from the analysis, select **Include**.
  - Green indicates that the loci are included. Gray indicates exclusion.
- 6. To modify the analysis thresholds, enter a value **0-100** in the applicable fields: Analytical Threshold, Interpretation Threshold, Stutter Filter, and Intralocus Balance.

The IT value must be greater than or equal to the AT value.

7. Select **Save** to apply the new settings.

## Delete an Analysis Method

You can delete any user-created analysis method. However, the default Verogen analysis methods cannot be deleted.

- 1. On the navigation rail, select **Settings**.
- Select Analysis Methods.
- 3. In the Select an Assay Type list, select the assay that corresponds to the analysis method you want to delete.
- 4. Select the applicable analysis method.
- 5. In the upper-right corner of the analysis method, select Delete Analysis Method.

Select OK to confirm the deletion.

#### Change the Default Analysis Method

- 1. On the navigation rail, select **Settings**.
- Select Analysis Methods.
- In the Select an Assay Type list, select the kit that corresponds to the analysis method you want to set as the default.
- 4. Select the applicable analysis method.
  - A green star identifies the current default.
- 5. Set the selected analysis method as the default for the assay type by selecting the empty star next to the analysis method name:

# **CODIS** Report Defaults

To expedite report creation, define default values to appear in the CMF header of every CODIS report. You can override these values on a per-report basis.

#### **Define CODIS Report Defaults**

- 1. On the navigation rail, select **Settings**.
- Select CODIS Report Defaults.
- 3. Select mtDNA CODIS Report for mtDNA samples or CODIS Report for STR samples.
- 4. In each of the following fields, enter a default value:
  - Destination Laboratory ORI

    —Identifier for the agency that processes the CODIS report.
  - Source Laboratory ORI—Identifier for the agency that processes the specimen.
  - Submission User ID—Identification of the person who generates the CODIS report.
- 5. Select **Save** to apply the default values.

## Set Specimen Categories

The Specimen Categories setting determines which specimen categories are available to assign to samples when creating a CODIS report. By default, all categories are available.

- 1. On the navigation rail, select **Settings**.
- Select CODIS Report Defaults.
- 3. Select mtDNA CODIS Report for mtDNA samples or CODIS Report for STR samples.
- 4. Select the Specimen Categories tab.
  - This tab lists all specimen categories, divided into CODIS and Custom columns.
- 5. Select a specimen category to make it available or unavailable for reports.
  - A blue toggle key indicates an available category.

#### Add Custom Specimen Categories

Add custom specimen categories to the software to supplement the predefined categories. Custom categories can be renamed, enabled or disabled, and deleted.

- On the navigation rail, select Settings.
- 2. Select CODIS Report Defaults.
- 3. Select mtDNA CODIS Report for mtDNA samples or CODIS Report for STR samples.
- 4. Select the Specimen Categories tab.
- 5. In the New Specimen Category Name field, enter a name of your choice.
- 6. Select Add.

The new category appears in the Custom column and is enabled by default.

# **User Management**

The Admin Settings page lists each account by user name. On this page, an administrator can manage permissions for each account, including granting access to new accounts and assigning access levels.

#### **Enable and Disable Accounts**

New accounts are disabled by default and must be enabled by an administrator.

- 1. Select Settings.
- 2. Select User Management.
- 3. Select Enabled to enable or disable an account.

Blue indicates an enabled account.

# Assign Access Levels

- 1. Select Settings.
- Select User Management.
- Select Administrator to enable or disable administrator access for an account.

Blue identifies an administrator account.

# System Settings

The System Settings page hosts the history setting. When enabled, this setting lets you view system events and user actions for runs, projects, and samples and includes these activities in project and sample reports.

By default, the history setting is enabled. However, the software *always records activity*, regardless of whether the setting is enabled. When enabling the setting, you can choose whether to store or delete activities that occurred since the last disabling.

#### Enable Run, Project, and Sample History

- 1. On the navigation rail, select **Settings**.
- 2. Select System Settings.
- 3. Select Enable or disable visibility of system events and user actions recorded by the software.
- 4. In the dialog box, select whether to store or delete previous system events and user actions:
  - **Permanently Delete All Previous History**—Start recording and delete all previous activity. You can access only activities recorded from this point forward.
  - Store Previous History—Continue recording and grant access to all activity starting from the creation of a run, project, or sample.

An Activity icon appears on the Run Details page, Project sidebar, and Results page. Selecting the icon opens the history.

## Disable Run, Project, and Sample History

- 1. On the navigation rail, select **Settings**.
- Select System Settings.
- 3. Select Enable or disable visibility of system events and user actions recorded by the software.

  The History icon disappears from the Run Details page, Project sidebar, and Results page.

## **Enable or Disable Dark Mode**

- 1. Select Settings.
- 2. Select **Dark Mode** to darken or brighten the user interface.

# Troubleshooting

Troubleshoot a Run	63
Troubleshoot Analysis Errors	66
Perform a Status Check	67

#### Troubleshoot a Run

The following sections provide recommendations for troubleshooting run problems. Most recommendations include library prep steps, which are documented in the library prep kit reference guides. Visit verogen.com/documentation to download the guide for your kit.

## Low Intensity

Low intensity (number of reads) can indicate a problem with sample quantity or quality.

- Make sure that DNA concentrations are not below the targeted amount.
- Depending on the source, purify or dilute the DNA and resequence the sample.
- Resequence with fewer samples, maintaining a minimum of eight samples per run.

#### **High Cluster Density**

Review run metrics: If cluster density is high, but sample representation shows that reads are above the sample read count guideline, proceed with analysis. The following factors can result in insufficient data for analysis.

- During normalization, library from the Purified Library Plate or the Normalization Storage (LNB1) were overpipetted.
  - If the sample has a high number of alleles with quality control indicators, be sure to use correct volumes and reagents for normalization and resequence the sample.
  - If the sample volume is correct, further dilute libraries in the Pooled Normalized Libraries tube.
  - If normalization volumes or reagents were incorrect, repeat normalization and resequence.
- Too much adapter dimer in library prep.
  - If your laboratory has a Bioanalyzer or Fragment Analyzer, analyze 1-2 μl library from the Purified Library Plate to check the amount of primer dimer.
    - If the checked library has > 5% adapter dimer, repeat library purification using the 30 μl library remaining in
      each well of the Purified Library Plate. Adjust purification to use 30 μl Sample Purification Beads 2 (SPB2),
      add 32.5 μl Resuspension Buffer (RSB), and transfer 30 μl to a new 96-well PCR plate. Proceed with
      normalization.
    - For questions, contact Verogen Technical Support.
  - [Optional] If SPB2 is limited, purify libraries in the Purified Library Plate using a standard column cleanup method. Follow the manufacturer instructions.
- An insufficient volume of Hybridization Buffer (HT1) was added, or too much normalized library was added to HT1.
  - If the sample has a high number of alleles with quality control indicators, resequence the sample. Make sure the normalized library is sufficiently diluted.

#### Low Cluster Density

Review run metrics: If cluster density for the run is high, but sample representation shows that reads are above the sample read count guideline, proceed with analysis. If the sample has an insufficient number of reads due to low density, the problem might be due to the following factors.

- DNA input was too low or too degraded, or the DNA dilution was incorrect. If possible, repeat library prep with more DNA.
- · A critical reagent was not added during library prep. Repeat the prep.
- LNB1 was insufficiently pipetted during normalization. Consider diluting libraries in the Pooled Normalized Libraries tube less.
- HT1 was overpipetted during library prep. Repeat sequencing using the correct volume.
- Excessive sample was lost during purification. Repeat the library prep using best practices for bead handling.
- After heating, the bead-based normalized library pool was not immediately added to HT1. If possible, repeat the prep.
- The library was not otherwise denatured properly.
  - Make sure that the final library pool is denatured with 2 N NaOH (HP3) that has a pH > 12.5.
  - Make sure that the final library pool is heat-denatured, if applicable.
- The MiSeq FGx System has a blockage.
  - Perform a maintenance wash and repeat sequencing.
  - Refill the wash tray and bottle for every wash.
  - Make all wash solutions daily.

# Low Percentage of Clusters Passing Filter

If the percentage of clusters passing filter is too low, the problem might be due to the following factors.

- The cluster density is too high. See High Cluster Density on page 63.
- The MiSeq FGx System has a blockage.
  - Perform a maintenance wash and repeat sequencing.
  - Refill the wash tray and bottle for every wash.
  - Make all wash solutions daily.
- If the problem remains unresolved, reagents might not be performing as expected. Contact Verogen Technical Support.

# High Phasing and Prephasing

The following factors can cause high phasing or prephasing.

- The cluster density is too high. See High Cluster Density on page 63.
- The ambient temperature is too high.

- Verify the appropriate temperature in the MiSeq FGx Sequencing System Site Prep Guide (document # VD2018012).
- Check the MiSeq FGx System temperature logs to confirm whether the ambient temperature is too high.
   Contact Verogen Technical Support for help locating the logs.
- The reagent chiller temperature is too low or too high. The system displays an error.
- The MiSeq FGx System has a blockage.
  - Perform a maintenance wash and repeat sequencing.
  - Refill the wash tray and bottle for every wash.
  - Make all wash solutions daily.
- If the problem remains unresolved, reagents might not be performing as expected. Contact Verogen Technical Support.

#### Incomplete Run

The following factors can cause a run to stop prematurely.

- The HSC was not added during library prep. Add HSC per denature and dilute instructions and repeat the run.
- Cluster density is too high. See High Cluster Density on page 63.
- · If an error message is displayed, contact Verogen Technical Support.

### Low Quality Scores for Read 1 or Read 2

The following factors can cause low quality scores for Read 1 and Read 2.

- Phasing and prephasing are too high. See High Phasing and Prephasing on page 64.
- The cluster density is too high. See High Cluster Density on page 63.
- If the problem remains unresolved, reagents might not be performing as expected. Contact Verogen Technical Support.

## Low Quality Scores for Index 1 Read or Index 2 Read

The following factors can cause low quality scores for the Index 1 Read or the Index 2 Read.

- Phasing and prephasing are too high. See Phasing and Prephasing on page 69.
- The cluster density is too high. See High Cluster Density on page 63.
- Not enough samples were sequenced, so diversity in the index reads was low.
- If the problem remains unresolved, reagents might not be performing as expected. Contact Verogen Technical Support.

## Low Reads per Sample

Review sample representation: If the samples exceed the sample read count (intensity) guideline, proceed with analysis. If the samples are below the guideline, the following factors are the likely cause.

DNA input was too low or the DNA dilution was incorrect. If possible, repeat library prep with more DNA.

- A critical reagent was not added during library prep. Repeat the prep.
- LNB1 was insufficiently pipetted during normalization.
  - If the sample has a high number of alleles with quality control indicators, resequence the sample.
  - If the sample has a low number of alleles with quality control indicators, consider diluting the libraries in the Pooled Normalized Libraries tube less and resequence the sample.
  - If the run does not pass, repeat library normalization.
- HT1 was overpipetted during library prep. Repeat library normalization.
- The libraries were not heat-denatured during library prep. Perform the heat-denature step in the Diluted Normalized Libraries tube and repeat the run.

## Low Reads per Sample for the HSC

Review sample representation: If the samples exceed the sample read count guideline, proceed with analysis. If the samples are below the guideline, the following factors are the likely cause:

- HSC was not added during library prep. Add HSC per denature and dilute instructions and repeat the run.
- HT1 was overpipetted during library prep. Repeat library normalization.
- The libraries were not heat-denatured during library prep. Perform the heat-denature step and repeat the run.
- HSC was not denatured with HP3. Make sure HSC is denatured with HP3 before adding it to the library, and then repeat the run.

# **Troubleshoot Analysis Errors**

The following sections provide recommendations for troubleshooting potential problems with the ForenSeq Analysis service or network.

# Analysis Error Message

When creating an analysis results in an error message, the application might not be running or the system cannot locate the run data for analysis. Use the following steps to troubleshoot.

- 1. From Windows Start on the server desktop, type Services, and then select **Services app**.
- 2. In the Services dialog box, locate ForenSeq Analysis.
- 3. Check the Status column:
  - If it is blank, right-click the service and select Start. After the service has started, retry analysis creation.
  - If it is running, complete the remaining steps.
- 4. Open the log file for the analysis at the following path:
  - C:\Illumina\Forenseq UAS\Analysis\logs\application.log
- 5. In the log file, find the message EXCEPTION OCCURRED: System. Exception.
  - This message confirms that the system cannot locate the run data.
- 6. If the run repository is stored in a network location, make sure that the network is running.

- 7. Navigate to the run folder listed in the log file as RunstoragePath to make sure that the data are in the expected location.
- 8. Retry analysis creation.

#### Quality Metrics Icons Are Unavailable

When all Quality Metrics icons on the Home page are unavailable, the system might be unable to locate data or the network where data are stored is disconnected. Use the following steps to troubleshoot.

- 1. From Windows Start on the server desktop, type Services, and then select **Services app**.
- 2. In the Services dialog box, locate ForenSeq Analysis.
- 3. Check the Status column:
  - If it is blank, right-click the service and select Start.
  - If it is running, complete the remaining steps to check network connectivity.
- 4. Navigate to the run folder described in the configuration file.
- 5. Make sure that data are in the expected location.
- 6. Make sure that no permission changes were made to access the run folder location.

#### Perform a Status Check

A status log in UAS v2.0 aggregates information from a run, sample, or installation failure and a history of actions. Use the following instructions to export this information to send to Verogen when troubleshooting.

- 1. On the navigation rail, select **Runs**, and then select a run.
- 2. On the Run Details page, select Run Details, and then select Run Status Check.
- In the Status Log Report dialog box, select Copy.
   The copy command copies the content displayed under System Status and Most Recent Log to your clipboard.
- 4. Paste the copied content to an external file.
- 5. When you are finished, select Close.

# **Run Metrics**

Accessing Run Metrics	68
Quality Metrics	68
Positive and Negative Controls	73
Sample Representation	75

# **Accessing Run Metrics**

Available in Projects after sequencing, run metrics provide run quality, positive and negative amplification control results, and sample representation for an overall view of run performance.


To open run metrics, select a QPN icon on the Project side sheet. A subset of these metrics, run quality and sample representation, are also available on Run Details pages.

# **Quality Metrics**

Quality metrics monitor run quality during sequencing. After sequencing, the software preserves the quality metrics and makes them available in both Runs and Projects. In Runs, the quality metrics appear in a pop-up window. In Projects, they appear on a page.

Color-coded bars indicate results for each metric and show optimum ranges for cluster density, clusters passing filter, phasing, and prephasing. Short, vertical bars with corresponding icons indicate the quality of each read and the HSC.

Figure 7 Quality metrics overview

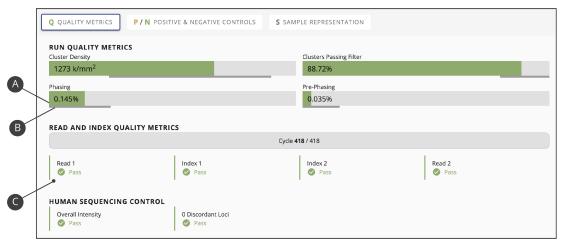

- A Color-coded horizontal bars and values indicate run quality metrics.
- B Thinner horizontal bars in dark gray indicate optimum ranges.
- C Color-coded vertical bars and icons indicate the quality of each read and the HSC.

## Run Quality Metrics

Run quality metrics include cluster density, clusters passing filter, phasing, and prephasing values. The following colors indicate overall quality.

Color	Indication
	Values are within the target range.
	Values require further evaluation of run quality.
	The target range for the metric.

#### **Cluster Density**

Cluster density is the number of clusters per square millimeter (K/mm²). For ForenSeq libraries, the target cluster density is 400-1650 K/mm². Values outside this range can still produce results sufficient for analysis. However, substantial deviations from the target range can impact other quality metrics and decrease the amount of usable data.

#### Clusters Passing Filter

Clusters passing filter is the percentage of clusters that passed the quality filter. This metric is based on the Illumina chastity filter, which measures quality and can detect low-quality base calls. Data appear after cycle 25.

The chastity of a base call is the ratio of the intensity of the greatest signal divided by the sum of the two greatest signals. If multiple base calls have a chastity value < 0.6 in the first 25 cycles, reads do not pass filter.

For ForenSeq libraries, the target value for a cluster passing filter is  $\geq$  80%. Values < 80% can still produce results that are sufficient for analysis. However, substantial deviations from the target can impact other quality metrics and decrease the amount of data.

#### Phasing and Prephasing

For both phasing and prephasing, lower percentages indicate higher quality run statistics. Phasing and prephasing values outside the target ranges can still produce results that are sufficient for analysis.

- Phasing shows the percentage of molecules in a cluster that fall behind the current cycle in Read 1. For ForenSeq libraries, the target phasing value is ≤ 0.25%.
- Prephasing shows the percentage of molecules in a cluster that jump ahead of the current cycle in Read 1. For ForenSeq libraries, the target prephasing value is ≤ 0.15%.

#### Read and Index Quality Metrics

The read and index quality metrics use the following colors and icons to indicate the status of each read and overall quality.

#### Table 9 Colors

Color	Indication
	The average quality of assessed reads is within the target range.
	The average quality of assessed reads is outside the target range.
	The read has not yet occurred.
Table 10 Icons	
Icon	Indication

Icon	Indication
	The average quality of assessed reads is within the target range.
•	The average quality of assessed reads is outside the target range.
	The read has not yet occurred.

#### Reads in a Run

A sequencing run completes up to four reads. Read 1 and Read 2 sequence the DNA template strands, and the Index 1 Read and Index 2 Read sequence the index adapters.

- Read 1—Read 1 sequencing primer is annealed to the template strand during cluster generation. RTA evaluates the first 50 cycles for quality.
- Index 1 Read—The Read 1 product is removed and the Index 1 sequencing primer is annealed to the same template strand as in Read 1. After index read preparation, the Index 1 Read is performed. RTA evaluates all eight cycles for quality.
- Index 2 Read—The Index 1 Read product is removed and the template anneals to the P5 primer grafted to the flow
  cell surface. The run proceeds through seven chemistry-only cycles without any imaging, followed by eight cycles of
  sequencing. RTA evaluates all eight cycles for quality.
- Read 2—The Index 2 Read product is extended to copy the original template strand. The original template strand is then removed and the Read 2 sequencing primer is annealed.

#### Cycles per Read

Quality metrics for each read appear after the read is complete. The number of cycles in each read depends on the assay.

Assay Type	Read	Cycles	Quality Metrics Appear
ForenSeq Kintelligence	Read 1	1-151	Cycle 152
	Index 1 Read	152-159	Cycle 160
	Index 2 Read	160-167	Cycle 168
	Read 2	168-318	After cycle 318
ForenSeq mtDNA Whole Genome	Read 1	1-201	Cycle 202
	Index 1 Read	202-209	Cycle 210
	Index 2 Read	210-217	Cycle 218
	Read 2	218-418	After cycle 418
ForenSeq mtDNA Control Region	Read 1	1-151	Cycle 152
	Index 1 Read	152-159	Cycle 160
	Index 2 Read	160-167	Cycle 168
	Read 2	168-318	After cycle 318
ForenSeq MainstAY	Read 1	1-351	Cycle 352
	Index 1 Read	352-359	Cycle 360
	Index 2 Read	360-367	Cycle 368
	Read 2	368-398	After cycle 398

## **Human Sequencing Control**

The software assesses the HSC and uses the following colors indicate the overall result.

Color	lor Indication	
	The control meets or exceeds the criteria.	
	The control does not meet at least one criterion.	

When the HSC is orange, evaluate the run quality metrics to help determine whether to repeat the run. If the run quality metrics are green, indicating that values are within target ranges, you can interpret the sample with caution.

#### Loci

The HSC includes the loci in the following Table. The software lists any loci that do not meet intensity or genotype concordance criteria.

Loci	151 or 201 Cycles	351 Cycles
D3S1358	<b>~</b>	<b>~</b>
D5S818	<b>~</b>	<b>~</b>
D7S820	~	<b>~</b>
D8S1179	<b>~</b>	<b>~</b>
D13S317	<b>~</b>	<b>~</b>
D16S539	<b>~</b>	<b>~</b>
D18S51	~	<b>~</b>
D21S11		<b>~</b>
FGA	<b>~</b>	<b>~</b>
PentaD	~	<b>~</b>
PentaE	~	<b>~</b>
TPOX	~	<b>~</b>
DYF387S1		<b>~</b>
DYS391	<b>~</b>	<b>~</b>
DYS392	~	<b>~</b>
DYS448		<b>~</b>
DYS460	<b>~</b>	<b>~</b>
Y-GATA-H4		<b>~</b>
DXS7423	<b>~</b>	<b>~</b>
DXS10074		<b>~</b>
DXS10103	<b>~</b>	<b>~</b>

# Positive and Negative Controls

The software assess the same positive and negative control metrics for every run, regardless of assay type. However, the assay type determines the type of positive amplification control and whether the software evaluates SNPs typed, STRs typed, or positions called.

Table 11 Control assessments

Assay Type	Positive Amplification Control	Value
ForenSeq Kintelligence	NA24385 Positive Amplification Control DNA (NA24385)	SNPs typed
ForenSeq mtDNA Control Region	Control DNA HL60 (HL60)	Positions called
ForenSeq mtDNA Whole Genome	Control DNA HL60 (HL60)	Positions called
ForenSeq MainstAY	NA24385 Positive Amplification Control DNA (NA24385)	STRs typed

### Positive Amplification Control

Positive amplification control metrics provide data for each sample identified as a ForenSeq Positive Amplification Control. Every sample is analyzed and compared to this control. The software then lists each sample with values for how many reads passed filter and SNPs typed, STRs typed, or positions called. Colors indicate overall success. If applicable, a table displays discordant positions or loci.

#### SNPs Typed

SNPs typed shows the number of loci with a base call out of the total number of loci for NA24385. A short, vertical bar displays the following colors to indicate call success and concordance with NA24385.

Color	Indication
	All loci have calls and are concordant.
	At least one locus does not have a call or contains a discordant call.

#### STRs Typed

STRs typed shows the number of loci with a base call out of the total number of loci for NA24385. A short, vertical bar displays the following colors to indicate call success and concordance with NA24385.

Color	Indication
	All loci have calls and are concordant.
	At least one locus does not have a call or contains a discordant call.

#### Positions Called

Positions called shows the number of positions, including indels, with a base call out of the total number of positions for HL60. A short, vertical bar displays the following colors to indicate call success and concordance with HL60.

	Color	Indication
All positions have calls and are concordant.		All positions have calls and are concordant.
		At least one position does not have a call or contains a discordant call.

#### Discordance Table

When a sample has calls that are not concordant with the positive amplification control calls, the positive control metrics include a discordance table. For each discordant call, the table displays the position number or locus where the discordance occurred, the call observed for the sample, and the expected call.

The table displays up to 50 positions or loci. When the number of discordant positions or loci exceeds 50, use the scroll arrows to move through the table.

#### Pass Filter Reads

Pass filter reads displays the following colors to indicate whether the number of reads that passed filter meet the guideline of ≥ 50,000 reads for ForenSeq mtDNA libraries, 15 million reads for ForenSeq Kintelligence libraries, or ≥ 15,000 reads for ForenSeq MainstAY libraries. Reads below the guideline might still be usable. Review quality metrics and sample data to help make a determination.

#### Table 12 Colors

Color	Indication
	The total number of pass filter reads exceeds the guideline.
	The total number of pass filter reads is below the guideline.
	The guideline for the number of reads that pass filter.

# Negative Amplification Control

Negative control metrics provide data for each sample identified as a negative amplification control or reagent blank. The software lists each of these samples with the number of SNPs typed, STRs typed, or positions called and a color-coded indicator of overall success.

#### SNPs Typed (Kintelligence)

SNPs typed shows how many loci a base call was made for. The short vertical bar uses the following colors to indicate call success.

Color	Indication	
	No loci have a base call.	
	At least one locus has a base call.	

## STRs Typed (MainstAY)

STRs typed shows how many loci a base call was made for. The short vertical bar uses the following colors to indicate call success.

Color	Indication
	No loci have a base call.
	At least one locus has a base call.

### Positions Called (mtDNA)

Positions called shows how many positions, including indels, a base call was made for. The short vertical bar uses the following colors to indicate call success.

Color	Indication
	No positions have a base call.
	At least one position has a base call.

# Sample Representation

Sample representation is a metric that provides quantitative sample and run information. After sequencing, the software preserves sample representation data and makes it available in both Runs and Projects. In Runs, sample representation data appear in a pop-up window. In Projects, these data appear on a page.

Sample representation shows an index CV percentage for the run and the number of reads (intensity) and read distribution for each sample. These data are presented in a bar chart that also shows the sample read count guideline. Hovering over a bar displays the exact number of reads for the sample.

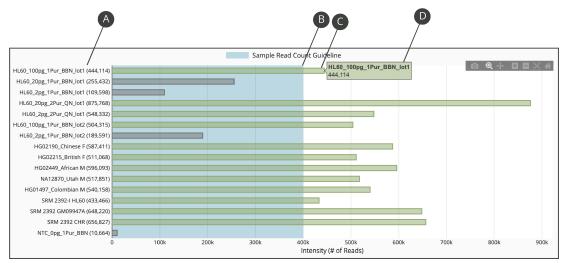



Figure 8 Sample Representation bar chart

- A Sample name and exact number of reads for the sample
- B Threshold of the read count guideline
- C Bar indicating number of reads for the sample
- D Exact number of reads for the sample (hovering)

#### Index CV

Above the sample representation bar chart, Index CV shows a percentage that represents the distribution of read counts for all samples in the run. CV is the coefficient of variation for the number of read counts across all indexes.

Figure 9 Index CV percentage

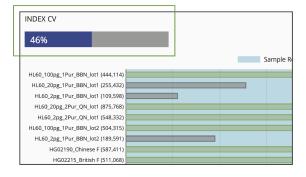

### Sample Read Count Guideline

The sample read count guideline varies by assay type. As guidelines, these values are intended to help with quality reviews of the run and samples. Samples with fewer reads can still demonstrate complete coverage and provide enough data for interpretation. Review quality metrics and sample data to help make a determination.

Assay Type	Sample Read Count Guideline
ForenSeq Kintelligence	15 million reads per sample
ForenSeq mtDNA Control Region	50,000 reads per sample
ForenSeq mtDNA Whole Genome	400,000 reads per sample
ForenSeq MainstAY	15,000 reads per sample

The following colors indicate how sample read counts compare to the applicable sample read count guideline.

Color	Indication
	The sample read count guideline.
	The sample meets or exceeds the guideline.
	The sample does not meet the guideline and requires careful interpretation.

# Report Formats

File Formats for mtDNA Reports	78
File Formats for Kintelligence Reports	80
File Formats for MainstAY Reports	83

# File Formats for mtDNA Reports

The following sections describe the format and contents of each type of report the software generates for mtDNA samples. For instructions on creating reports, see *Report Generation* on page 39.

#### **CODIS** Report

A CODIS report is a .cmf file that includes a CMF header, reference sequence, and specimen profile. The specimen profile provides the following information for each sample in the report:


- Specimen ID
- Specimen category
- · Source ID
- Fragments

#### **EMPOP Report**

An EMPOP report is a .txt file. For each sample, the report includes a sequence string in FASTA format, variant calls in rCRS format, and names for the following items:

- Project
- Sample
- Assay
- · Analysis method
- Run

Figure 10 Example EMPOP report

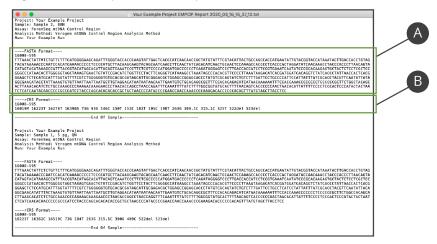

- A Sequence string in FASTA format
- B List of variant calls in rCRS format

### **Project Report**

An mtDNA project report is an Excel workbook. The first worksheet, Variant Sample Report, provides the following project information:

- Project name and optional description
- Project creation date and time
- · Run and sample counts
- Sample data:
  - Single nucleotide variant (SNV) count per sample
  - Insertion count per sample
  - Deletion count per sample
  - Call or no call count per sample
  - Manually changed calls per sample
  - Variants per sample

If applicable, the second worksheet in the project report, Project History, includes administrative details for the project and a project history log that lists time-stamped system events and user actions.

Each subsequent worksheet contains sample information (one worksheet per sample). For details, see Sample Report.

## Sample Report

An mtDNA sample report is an Excel workbook. If applicable, the first worksheet is Sample History, which includes administrative details for the sample and a sample history log that lists time-stamped system events and user actions.

The second worksheet is named for the sample and provides the following information:

- Assay type
- · Analysis method
- Run name
- Analytical and interpretation thresholds
- Minimum Q-score
- Minimum read count
- SNVs
- · Insertion and deletion counts
- Call count for reagent blanks and negative controls
- No call count for samples and positive controls
- Manually changed calls
- Variant details:
  - Variant and reference calls
  - Total read count
  - Strand depth
  - Reads A, C, T, and G
  - Percent A, C, T, and G
  - Reads del
  - Percent del
  - Reads ref
  - Percent ref

## File Formats for Kintelligence Reports

The following sections describe the format and contents of each type of report the software generates for Kintelligence samples. In the reports, rsID numbers identify individual SNPs. For instructions on creating Kintelligence reports, see *Generating Kintelligence Reports* on page 43.

## Phenotype and Ancestry Report

A phenotype and ancestry report is an Excel workbook comprised of four worksheets that detail the estimates of hair color, eye color, and biogeographical ancestry: Estimation, SNP Data, SNP Balance Figure, and Settings. In the Estimation and SNP Data worksheets, a header section lists the software version, analysis method, sample gender, project creation date and time, the user who generated the report, and sample, project, and run names.

- Estimation
  - Hair color results: brown, red, black, and blonde
  - Eye color results: intermediate, brown, blue
  - Biogeographical ancestry results (plot): the sample with distance to the nearest centroid

- Population information for samples in the centroid: 1000 Genomes population and abbreviation, count, and data used to train models
- SNP Data
  - Allele name, whether the locus is typed, and the number of reads for each locus
  - Number of loci typed, genotype, and QC indicators for hair and eye color SNPs, common SNPs, and biogeographical ancestry SNPs

Loci typed is the number of aSNP and pSNP loci identified out of total targeted aSNP and pSNP loci.

- SNP Balance Figure
  - Whether the allele is typed
  - Number of reads per base for each SNP with a grand total
- Settings
  - Software version
  - AT, IT, and intralocus balance for each locus

#### **Estimating Phenotype**

The software uses HIrisPlex, a multinomial logistic regression model, to determine the probabilities for hair and eye color. The software publishes these probabilities to the Estimation worksheet. 1-3

<sup>1</sup> Walsh, Susan, Lakshmi Chaitanya, Lindy Clarisse, Laura Wirken, Jolanta Draus-Barini, Leda Kovatsi, Hitoshi Maeda, et al., "Developmental validation of the HIrisPlex system: DNA-based eye and hair colour prediction for forensic and anthropological usage," *Forensic Science International: Genetics* 9 (March 2014): 150-161, https://doi.org/10.1016/j.fsigen.2013.12.006.

<sup>2</sup> Liu, Fan, Katevan Duijn, Johannes R. Vingerling, Albert Hofman, André G. Uitterlinden, A. Cecile J.W. Janssens, and Manfred Kayser, "Eye color and the prediction of complex phenotypes from genotypes," *Current Biology* 10, no. 5 (March 2009): R192-R193, https://doi.org/10.1016/j.cub.2009.01.027.

<sup>3</sup> Walsh, Susan, Fan Liu, Kaye N. Ballantyne, Mannis van Oven, Oscar Lao, and Manfred Kayser, "IrisPlex: A sensitive DNA tool for accurate prediction of blue and brown eye colour in the absence of ancestry information," *Forensic Science International: Genetics* 5, no. 3 (June 2011): 170-180, https://doi.org/10.1016/j.fsigen.2010.02.004.

## **Estimating Biogeographical Ancestry**

The software uses principal component analysis (PCA) to determine estimates of biogeographical ancestry. Data from Phase I of the 1000 Genomes Project (1000genomes.org) was used to train the model on the following super populations: European, East Asian, and African (excepting African Ancestry in Southwest USA [ASW]). For context, the unknown sample is projected with the Ad-Mixed Americans super population onto the first two pretrained components based on aSNP genotype calls.<sup>4</sup>

On the Estimation worksheet, the distance to the nearest centroid indicates how related a sample is to the general grouping for the centroid. For comparison, the worksheet displays the distance for 1000 Genome Project samples that contribute to the centroid. The chart on the Estimation worksheet includes one centroid for each of the major ancestries and orientational centroids at the one-quarter intervals between these groups. Centroids contextualize results and logically group populations.

<sup>4</sup> Abecasis, Gonçalo R., David Altshuler, Adam Auton, Lisa D. Brooks, Richard M. Durbin, Richard A. Gibbs, Matt E. Hurles, et al., "A map of human genome variation from population-scale sequencing," *Nature* 467 (October 2010): 1061-1073, https://doi.org/10.1038/nature09534.

### Generating SNP Data

The SNP Data worksheet provides data for aSNPs and pSNPs. Estimating hair and eye color requires genotype results for hair color, eye color, and common SNPs. Common SNPs are SNPs that estimate hair color, eye color, and biogeographical ancestry.

If one of these SNPs is not typed, the software does not generate a result. In contrast, generating a result for biogeographical ancestry requires that only one biogeographical ancestry SNP be typed. The absence of SNPs in the biogeographical ancestry estimation diminishes the accuracy of the estimation. With each execution, the estimation model is retrained using only the aSNPs with a multilocus genotype.

#### **GEDmatch PRO Report**

A GEDmatch PRO report is a .txt file. Using human genome assembly GRCh37 (hg19) as the reference genome, the report lists each SNP locus call by rsID with the following information:

- Chromosome and position
- Allele 1 and allele 2 typed calls reported on the plus strand

## **Project Report**

A Kintelligence project report is an Excel workbook. A header section lists the project name, project creation date and time, and the user who generated the report. A Samples table lists each sample in the report with the following information:

- · Analysis method
- Contributor status
- Gender
- Locus call rate
- Sample name and optional description
- · Total sample reads

#### Sample Report

A Kintelligence sample report is an Excel workbook with one worksheet per SNP type: X-, Y-, ancestry, phenotype, identity, and kinship. The workbook also includes a Settings worksheet and, if applicable, a Sample History worksheet. The Settings worksheet lists AT, IT, and intralocus balance for each locus. The Sample History worksheet includes administrative details for the sample and a sample history log listing time-stamped system events and user actions.

The header section of each SNP worksheet lists the analysis method, gender, contributor status, project creation date and time, user who generated the report, and the sample, project, and run names. Coverage information appears under the header, followed by a table listing information for each locus:

- Locus
- Genotype with allele name
- Whether the allele is typed
- Reads
- Quality control indicators

## File Formats for MainstAY Reports

The following sections describe the format and contents of each type of report the software generates for MainstAY samples. For instructions on creating reports, see *Report Generation* on page 39.

## **CODIS** Report

A CODIS report is a .cmf file that includes a CMF header, reference sequence, and specimen profile. The specimen profile provides the following information for each sample in the report:

- Specimen ID
- Specimen category
- Source ID
- Fragments

## Project Report

A MainstAY project report is an Excel workbook with two worksheets each for aSTRs and Y-STRs. All four worksheets include a header section that lists the project name, project creation date and time, and user who generated the report.

- The Autosomal STRs and Y STRs worksheets provide genotype information formatted in a table that lists the following information for each sample in the project.
  - Sample name and optional description
  - Analysis method
  - Contributor status
  - Gender
  - Typed alleles for each locus in the sample

- The Autosomal STR Coverage and Y STR Coverage worksheets provide coverage information formatted in a table that lists the following information for each sample in the project.
  - Sample name
  - Analysis method
  - Locus name and each allele name
  - Whether each allele is typed
  - Number of reads for each allele
  - Repeat sequence of each allele

### Sample Report

A MainstAY sample report is an Excel workbook with two worksheets each for aSTRs and Y-STRs and a settings worksheet that lists analysis thresholds for each locus.

- The Autosomal STRs and Y STRs worksheets include the following sample information.
  - A header section lists the following information for the sample:
    - Sample name
    - Project name
    - · Analysis method
    - Run name
    - Gender
    - Contributor status
    - Project creation date and time
    - User who generated the report
  - Locus information is formatted in a table that lists each locus with genotypes and any QC indicators. A fraction
    indicates how many loci are typed out of the total possible number of loci.
  - Coverage information is formatted in a table that lists the following information for each locus:
    - · Locus name and each allele name
    - Whether each allele is typed
    - Number of reads for each allele
    - Repeat sequence of each allele
- The Autosomal STR Figure and Y STR Figure worksheets indicate the total number of reads for the sample with a breakdown of how many reads each allele received.
  - This information is presented in both a table and a bar chart.
  - Filters for locus name and typed or untyped alleles are available at the top of the table.
- The Settings worksheet displays the analysis module version and lists the AT, IT, stutter filter, and intralocus balance for each locus in the sample.

# **Supporting Information**

Special Handling of mtDNA Data	85
Contributor Status of Kintelligence Samples	89
Gender Estimation for Kintelligence Samples	89
Contributor Status for MainstAY Samples	89
Gender Estimation for MainstAY Samples	89
Stutter filtering for MainstAY Samples	90

# Special Handling of mtDNA Data

The following sections describe how the software handles special circumstances, such as ambiguous base call codes, when analyzing sequencing data from mtDNA samples.

#### Consensus Calls

The software supports variant calling for multiple bases at a DNA position. When multiple bases exceed the interpretation threshold, the software reports the International Union of Pure and Applied Chemistry (IUPAC) code as the consensus for the position.

IUPAC Code	Base
А	Adenine
С	Cytosine
G	Guanine
Т	Thymine
R	A or G
Υ	C or T
S	G or C
W	A or T
K	G or T
M	A or C
В	C, G, or T
D	A, G, or T
Н	A, C, or T
V	A, C, or G

## **Ambiguous Calls**

When data for the indel and a base call at a position exceed the IT, the software reports the ambiguous call code for an indel as the consensus for the position.

Base	Call
а	Adenine
С	Cytosine
g	Guanine
t	Thymine
r	a or g
у	cort
S	g or c
W	a or t
k	gort
m	a or c
b	c or g or t
d	a or g or t
h	a or c or t
V	a or c or g
n	Any base

## Ambiguous Call Example

An example of an ambiguous call, the following figure shows an ambiguous deletion at position 501. Because an equal number of reads support a deletion call and a reference call of G, the software reports an ambiguous call of g as the consensus.

Figure 11 Ambiguous deletion at position 501

Coordinate	498	499	500	501	502	503
Call	Α	С	Т	g	Α	С
Reads	A	С	Т		Α	С
	A	С	Т		Α	С
	A	С	Т	_	Α	С
	A	С	Т	_	Α	С
	A	С	Т		A	С
	A	С	T		A	С
	A	С	Т		Α	С
	A	С	T		Α	С
	A	С	T	_	Α	С
	A	С	Т	G	Α	С
	A	С	Т	G	Α	С
	A	С	Т	G	Α	С
	A	С	Т	G	A	С
	A	С	Т	G	A	С
	A	С	Т	G	A	С

#### **SWGDAM Nomenclature**

The software follows SWGDAM nomenclature. For variant calling, differences between the reference and sample sequences are recorded as polymorphisms with the nucleotide position and DNA base difference indicated in the haplotype nomenclature (for example, 16182C).

### **Primer Trimming**

The software trims primer sequences for ForenSeq mtDNA amplicons, which are based on a tiled primer design. This trimming allows unambiguous interrogation of mtDNA amplicon inserts without degenerate oligonucleotides. The mtDNA Variant Processor v1.0.0 application in Illumina BaseSpace Sequence Hub does not trim primers from amplicons generated from a tiled primer design.

#### Nuclear mtDNA Insertions

Nuclear DNA sequences homologous to mtDNA are often called nuclear mtDNA insertions (NUMTs). During library prep, these sections can be enriched alongside the respective mtDNA target. When performing alignment, UAS v2.0 identifies reads that map closely to sequences in the Reference Human NumtS (RHNumtS) compilation and removes these reads before further analysis. <sup>5,6</sup>

<sup>&</sup>lt;sup>5</sup> Mishmar, Dan, Eduardo Ruiz-Pesini, Martin Brandon, and Douglas C. Wallace, "Mitochondrial DNA-like Sequences in the Nucleus (NUMTs): Insights Into Our African Origins and the Mechanism of Foreign DNA Integration," *Human Mutation* 23 no. 2 (February 2004): 125-133, https://doi.org/10.1002/humu.10304.

<sup>&</sup>lt;sup>6</sup> Calabrese, Francesco Maria, Domenico Simone, Marcella Attimonelli "Primates and Mouse NumtS in the UCSC Genome Browser." Supplement, *BMC Bioinformatics* 13, no. S4 (March 2012): S15. https://doi.org/10.1186/1471-2105-13-S4-S15.

#### **HVII C-Stretch Coverage**

The C-stretch describes homopolymers of cytosine-stretches found in hypervariable region I (HVI) and hypervariable region II (HVII) of mtDNA. Primer pairs that amplify the region between position 262 and position 353 achieve coverage of the C-stretch (positions 303-315) in HVII. Sequences generated from the forward template strand have high accuracy and alignment. Sequences generated from the reverse template strand have high accuracy and alignment until reaching the C-stretch. Therefore, reads from the reverse strand that start sequencing at position 262 and do not meet alignment requirements are soft-clipped in UAS v2.0 after position 303. The soft-clipped positions are not used for base calling. As a result, approximately half the coverage (read counts) are obtained for positions 304-353 relative to positions 262-303 (on the reverse strand).

In the following figure, a Results page shows similar total read count and strand depth at positions 262-303 and 304-353. Positions 262-303, where reads from both strands are reported, have similar coverage. Total reads for positions 304-353 are about half that of positions 262-353.

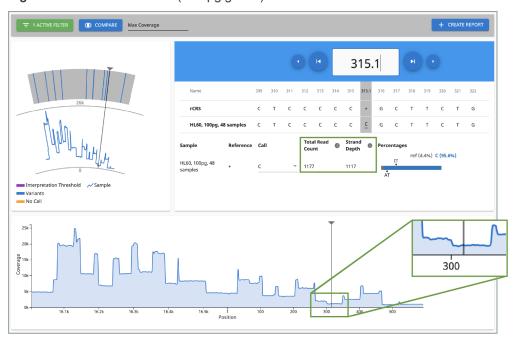


Figure 12 Results for HL60 (100 pg gDNA) in UAS v2.0 v2.1

<sup>&</sup>lt;sup>7</sup> Nakamura, Kensuke, Taku Oshima, Takuya Morimoto, Shun Ikeda, Hirofumi Yoshikawa, Yuh Shiwa, Shu Ishikawa, et al., "Sequence-specific error profile of Illumina sequencers," *Nucleic Acids Research* 39, no. 13 (July 2011): e90. https://doi.org/10.1093/nar/gkr344.

<sup>&</sup>lt;sup>8</sup> Allhoff, Manuel, Alexander Schönhuth, Marcel Martin, Ivan G. Costa, Sven Rahmann, and Tobias Marschall. Discovering motifs that induce sequencing errors," Supplement, *BMC Bioinformatics* 14, no. S5 (April 2013): S1. https://doi.org/10.1186/1471-2105-14-S5-S1.

# Contributor Status of Kintelligence Samples

The software applies the following logic to determine the contributor status of a Kintelligence sample.

- Any of the following scenarios indicate a mixture:
  - The sample has ≥ 10 Y-SNPs with two alleles.
  - The sample has ≥ 10 Y-SNP loci typed and ≥ 10 X-SNPs with two alleles.
  - The sample has zero Y-SNP loci typed with ≥ 60% autosomal SNP (auSNP) heterozygosity and ≥ 50%
     X-SNP heterozygosity.
- If the sample has ≤ 50% auSNP heterozygosity, ≤ 40% X-SNP heterozygosity, ≥ 65% auSNP intralocus balance, and zero Y-SNPs with two alleles, the sample is single-source.
- If the sample is a negative amplification control, the call rate is < 50%, or the scenarios for mixture and single-source do not apply, the contributor status is inconclusive.

To adjust heterozygosity percentages for call rate, the observed heterozygosity is divided by the sample call rate.

## Gender Estimation for Kintelligence Samples

The software applies the following logic to estimate the gender of a Kintelligence sample.

- If the sample has ≥ 10 Y-SNP loci typed, the gender is male (XY).
- If the sample has zero Y-SNP loci typed and the call rate is ≥ 50%, the gender is female (XX).
- If the sample is a negative amplification control, the contributor status is mixture, or the scenarios for male and female do not apply, the gender is inconclusive.

## Contributor Status for MainstAY Samples

The software applies the following logic to determine the contributor status of a sample processed with ForenSeq MainstAY.

- If the sample is not a negative amplification control, and no loci have more alleles called than possible for that locus the contributor status is single-source.
- If the sample is not a negative amplification control and atleast 3 loci have more alleles than possible for the locus the contributor status is a mixture.
- If the sample is a negative amplification control or 1 -2 of the loci have more alleles called than is possible for that locus, the contributor status is inconclusive.

# Gender Estimation for MainstAY Samples

The software applies the following logic to estimate the gender of a sample processed with the ForenSeq MainstAY kit. Gender is estimated prior to estimating determining Contributor Status.

• If the sample is not a negative control, atleast 50% of the autosomal STRs have signal above the analytical threshold and no Y-STR loci have signal greater than the analytical threshold, the gender is female (XX).

- If the sample is not a negative control and atleast 5 Y-STR loci have signal above the analytical threshold, the gender is male (XY).
- If the sample is a negative amplification control, the contributor status is a mixture or if 1-4 Y-STRs have an allele above the analytical threshold the gender is inconclusive. In additiona, if less than 50% of the autosomal STRs have signal above the analytical threshold and no Y-STR loci are typed the gender is inconclusive.

## Stutter filtering for MainstAY Samples

A sequence is considered stutter, with offset (n) if the sequences differ in length by k times the reference repeat unit length, the number of repeat units in the potential stutter are identical to an adjacent unit of length of the reference repeat unit and there are no other mismatches between the sequences. The software assesses differences in the sequences for k = 2, -1 and +1 repeat units. The stutter filter setting, which depends on the differences in repeat units is applied as is for the n-4 and n-5 positions. The stutter filter setting is squared at the n-8 or n-10 and n+4 or n+5 positions. The n+1 repeat stutter filter is nor adjustable by users, and is set as the square of the n-1 repeat stutter filter, with the exception of DYS460 where the n-1 filter is set as the square root of the n-1 repeat stutter. The expected stutter intensity is the product of the stutter filter setting and the called allele intensity at a repeat length. The stutter quality control indicator for elevated stutter is displayed under the following conditions.

- Uncalled read intensity, with sequence of potential stutter of a called allele.
- Uncalled read intensity is greater than the maximum expected stutter (%) if the potential parent allele.

# Index

1 1000 Genomes Project 80-81 A	base calls 14, 74-75, 85-86 modified 32 multiple 85 original 32	coverage plot 13, 31 filters 31 custom assays 26 customer support 95
activity setting 12, 61 adapter dimer 63 administrators access 61 accounts 9 privileges 55, 61 settings 61 alignment 87-88 allele calls 55, 82 alleles 16, 35, 63, 65, 82 modifying 16, 18 typing 36 untyping 36 ambiguous calls 85 amplicons 87 lengths 18 modifying 56-57 analysis files 50 analysis methods 26, 37 default settings 55 editing 56 favorites 60 filtering 58-59 sorting 58-59 viewing results 38 analysis parameters 55, 76 analysis status 19 analytical thresholds 14, 56 ancestry estimates 81 ancestry informative SNPs 16 ancestry reports 39, 43, 45 applications 8, 87 assay types 26, 55-56 AT 57-58, 80, 83-84	base pairs 18 BaseSpace Sequence Hub 87 biological sex 36, 89-90  C C-stretch 88 Call and No Call filters 32 call rates 89-90 calls 13, 16, 18, 31 cards 16, 33 centroids 81 chastity filters 69 Chrome 8 chromosomes 36, 82 cluster density 63-65, 68-69 clusters passing filter 64, 68-69 CMF headers 39-40, 45-46, 60, 78, 83 CODIS reports 39, 43, 45 color coding 10, 16-18 colors 29-30, 32 commenting 35, 51, 54 comments deleting 54 editing 54 comparison samples 32 compatible products 7 concordance 74 configuration, software 55 contact information, Verogen 95 contributor status 89-90 contributors 16, 18, 35, 39-40, 45-46, 83-84 control software 8 controls 26, 30, 68 comparing 49 metrics 73	cycle numbers 70  D  dark mode 55 dashboard 9-10 dashes, position viewer 14 data recording 62 database reports 39, 43, 45 default analysis settings 56 default views 33 deletions 86 differences, samples 32 discordance 30, 74 diversity 65 DNA inputs 63-65 DNS 8 documentation 95 Domain Name System 8 downloading reports 13, 39-48  E  EMPOP file format 41 EMPOP reports 39, 43, 45 enrichment 87 error messages 65 estimating ancestry 81 estimating eye color 81 estimating hair color 81 Excel workbooks 79-80, 82-84 exceptions 85  F  FASTA file format 78 female samples 89-90 file formats 80 filtering analysis methods 58-59 reports 13
BAM files 50 base call codes 85-86	types 73 coordinates 14 coverage 13, 83-84	samples 31, 33 filters 13, 18, 27, 32 Call and No Call 32

flow cells 7, 25 compatibility 7	internet requirements 8 InterOp folder 30	negative amplification controls 89-90 negative amplification controls 13,
editing type 25	interpretation thresholds 13-14, 56,	74
ForenSeq Analysis service 66	85-86	networks
G	intralocus balance 16, 18, 80, 83-84	requirements 8
GEDmatch PRO reports 39, 43, 45	IT 57-58, 80, 83-84	troubleshooting 66
gender 16, 18, 35	IUPAC codes 85	no calls 13-14, 31
genotype concordance 71	K	non-directed sequencing 56
genotypes 83	K/mm <sup>2</sup> 69	nucleotides 85
H	kinship SNPs 16	number of reads 15, 63 NUMTs 87
	1	
heatmaps 18	L	O
help, technical 95	library prep kits 7, 26	original calls 32
heterozygosity 89	reference guides 63	outputs 78, 83
HIrisPlex 81	library types 56-57	P
History icon 62	loci 16, 71	•
history setting 11-12, 51, 61	editing 57-58	pass filter reads 30, 75-76
home-brew assays 26, 56	searching 58-59	passwords, lost 9
home page 9	typed 16, 80, 82	paused runs 19
hotspots 14 HSC 65-66	user actions 54	permissions 61
loci 71	locus call rates 16, 18	phasing and prephasing 30, 65, 68-
quality metrics 29, 68	locus calls 16, 18 locus details 16, 36	69
results 71	logging in 9	phenotype and ancestry reports 39,
hyperlink run names 9	logging out 9	43, 45
hyperlink sample names 9	lost passwords 9	phenotype informative SNPs 16
hypervariable regions 88	•	phenotype reports 39, 43, 45 plexities 7, 65
I	M	populations 80
I	male samples 89-90	position viewer 13, 32
i5 read 70	matching samples 32	dashes 14
i7 read 70	maximum samples 7	filters 31
icons 10, 12, 17-19, 69	metrics, colors 29-30	positions 13, 15, 82
identical samples 32	minimum quality score 56-57	selecting 14, 39, 45
identity informative SNPs 16	minimum read count 56-57	positions called 30
importing samples 26	MiSeq FGx Control Software 8	positive amplification controls 13, 73
inconclusive contributor status 89	MiSeq FGx reagents 7	primer-directed sequencing 56
inconclusive gender 89-90	MiSeq FGx System 7-8, 19	primer pairs 88
indels 86	mixtures 16, 18, 36, 89-90	primers 70
Index 1 and Index 2 Reads,	models 81	Print icon 50
quality 69	mtDNA navigator 13, 31	printing 12
index adapters 26, 70	control region 13	privileges, administrator 61
index CV 75 Information icon 50	filters 31 multiple bases called 85	project reports 39, 43, 45
installing 9		project summary, printing 50
intensities 18	N	projects 10
intensity 16, 18, 65, 69, 71, 75	navigation rails 10	assigning samples 26, 28-29 creating 25-26, 29

descriptions 49 icons 12 naming 12, 49 summaries 9  Q Q icon 10 QC indicators 13, 16-18, 33 QPN 12, 68 QPN icons 13 quality control 63, 76 quality filters 69 quality metrics 10-11, 13, 15, 29, 70  R rails 10 rCRS 14, 78 read counts 14, 31, 76 read cycles 70 read definitions 70 read distributions 75 read numbers 31 Reads 1 and 2, quality 69 reagent blanks 74 reagent kits 7, 25 reagent performance 64-66 Real-Time Analysis 8 reanalyzing samples 19, 38 recording data 62 reference genomes 82 reference sequences 78, 83 reports contents 78, 80, 83 downloading 13, 39-48 filtering 13 finding 13 types 39, 43, 45 worksheets 79-80, 82-84 resolution 8 results 13, 16-17	run details 9 run metrics 64, 71 run monitoring 11, 29, 68, 75 run status 19 runs 10 general information 11 metrics 63, 71 performance 11, 29 quality 30, 68 status 29  S safety data sheets 95 sample comparisons, reporting 49 sample information 13, 16-17, 25, 39, 45, 79-80, 82-84 Sample Overview section, printing 50 sample position, changing 13 sample read count guideline 30, 63- 66, 75 sample reports 39, 43, 45 sample representation 11, 13, 29- 30, 63-66, 68 sample results 13, 16-17 sample sheets 25-26 entering UDIs 26 samples adding information 25 analyzing 13 assigning to projects 25-26, 29, 37, 49 differences 32 filtering 31, 33 importing 26 intensity 63 locating 27 modifying 28-29 no differences 32	Samples page printing 50 scatterplots 18 screen resolution 8 SDS 95 searching loci 58-59 reports 13 samples 27 selected comparison samples 32 sequences 14, 78 sequencing processes 70 servers 7, 50 sex chromosomes 36 single-source samples 16, 18, 36, 89 SNP types 16 SNPs calls 82 types 16 soft clipping 88 software reset 9 sorting 13, 27, 33, 58-59 specimens 39-40, 45-46, 78, 83 star icons 60 stopped runs 19 strand depth 14-15 stutter 18, 84 stutter filters 56, 59 super populations 81 support, library prep 63 system dashboard 9-10 system events, definition 61 T technical assistance 30, 95 technical notes 55 third-party assays 26, 56 thresholds 55-58, 76, 80, 83-84
results 13, 16-17		thresholds 55-58, 76, 80, 83-84 tiling 87
filtering 31, 33	removing from runs 29	total read count 15
Results page 9	searching 27	total reads 14, 16, 18
printing 50	_	troubleshooting files 30
rs numbers 82	sorting 33	typing alleles 16, 18, 36
rsID 82	status 19	typing loci 16
RTA 70	samples of interest 32	typing loci To
run cycles 70		

U

UDIs
definition 24
sample sheets 26
uninstalling 9
untyping alleles 16, 18
user access 61
user actions, definition 61



validations 56 variant calls 55, 85 variants 13-14, 78 VCF files 50

## W

washes 64 web browser 8

# X

XX gender 89-90 XY gender 89-90

# **Technical Support**

For technical assistance, contact Verogen Technical Support.

	Contact Information
Address	11111 Flintkote Avenue San Diego CA 92121 USA
Website	www.verogen.com
Email	techsupport@verogen.com
Telephone	+1.833.837.6436 toll-free (North America) +1.858.285.4101 tel +44 (0) 208 054 8706 (United Kingdom)

Safety data sheets (SDS)—Available for download from verogen.com/product-documentation.

 $\textbf{Product documentation} - \textbf{A} vailable for download from verogen.com/product-documentation.}$ 



## Meet any challenge

Verogen is a dedicated developer of human identification products for sequencing and analysis of forensic genomic samples. Working closely with the forensics community, Verogen places exceptional value on flexible, scalable solutions that deliver results when you need them most.

Verogen

+1.833.837.6436 toll-free (North America)

+1.858.285.4101 (outside North America)

techsupport@verogen.com

www.verogen.com