
January 2022

EZ1&2™ ccfDNA Handbook

For automated purification of circulating cell-free DNA (ccfDNA) from plasma or serum using the EZ1® Advanced XL or the EZ2 Connect instruments

Contents

Kit Contents	4
Shipping and Storage	5
Intended Use	5
Safety Information.....	6
Quality Control.....	6
Introduction.....	7
Principle and procedure	7
Adsorption to anion exchange magnetic particles	8
Washing of bound nucleic acids	8
Elution of purified nucleic acids.....	8
Automation.....	9
Starting material.....	9
Sample volumes	10
Sample numbers.....	10
Yield and size of nucleic acids.....	10
Equipment and Reagents to Be Supplied by User	12
Important Notes.....	13
Protocol: Purification of ccfDNA from up to 8 ml Serum or Plasma Using the EZ1&2 ccfDNA Kit.....	21
Troubleshooting Guide	25
References	27
Appendix A: Recommendations for Plasma Separation and Storage.....	28
Appendix B: Example of an EZ1 Advanced Report File	30

Ordering Information	33
Document Revision History	35

Kit Contents

EZ1&2 ccfDNA Kit	(48)
Catalog no.	954854
Number of preps	48
Large Volume Tube (7 ml)	4 x 24
Magnetic Bead Suspension EZ	48 x 0.075 ml
Elution Buffer EZE*	4 x 15 ml
Elution Tubes 1.5 ml	50
Reagent Cartridges, EZ1&2 ccfDNA†	48
Disposable Tip Holders	50
Disposable Filter-Tips	50
Q-Card‡	1
Quick-Start Protocol§	1

* Elution buffer EZE is supplied as a concentrate – the working solution is prepared by the EZ1 or EZ2 instrument during the automated procedure.

† Contains sodium azide as a preservative. See page 6 for safety information.

‡ The information encoded in the bar code on the Q-Card is needed for reagent data tracking using EZ1 Advanced XL instruments.

§ To view this handbook's corresponding Quick-Start Protocol, go to www.qiagen.com/HB-2891.

Shipping and Storage

EZ1&2 ccfDNA Kits are shipped at ambient temperature. Upon receipt, store all kit components dry at room temperature (15–25°C). Do not freeze the reagent cartridges. When stored properly, buffers and reagent cartridges are stable until the expiration date printed on the kit box lid.

Partially used bottles of elution buffer EZE can be stored for a maximum of 4 weeks.

Intended Use

EZ1&2 ccfDNA Kits are intended for molecular biology applications. These products are not intended for the diagnosis, prevention, or treatment of a disease.

All due care and attention should be exercised in the handling of these products. We recommend all users of QIAGEN® products to adhere to the NIH guidelines that have been developed for recombinant DNA experiments, or to other applicable guidelines.

Safety Information

When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, please consult the appropriate safety data sheets (SDSs). These are available online in convenient and compact PDF format at www.qiagen.com/safety where you can find, view, and print the SDS for each QIAGEN kit and kit component.

If liquid containing potentially infectious agents is spilt on the EZ1 Advanced XL or EZ2 Connect, please refer to the instrument user manual for decontamination instructions.

Quality Control

In accordance with QIAGEN's ISO-certified Quality Management System, each lot of the EZ1&2 ccfDNA Kits is tested against predetermined specifications to ensure consistent product quality.

Introduction

Circulating, cell-free, DNA (ccfDNA), such as tumor-specific extracellular DNA fragments or fetal DNA in maternal blood, is present in serum or plasma usually as short fragments, <1000 bp. The concentration of cfDNA in biological fluids such as plasma or urine is usually low and varies considerably between individuals, usually from 1–100 ng/ml.

EZ1&2 ccfDNA Kits enable efficient purification of these cfDNA fragments from human plasma or serum samples. These samples can be either fresh or frozen (provided that they have not been frozen and thawed more than once). Human serum or plasma samples can be generated using blood collection tubes such as the PAXgene® Blood ccfDNA Tube. cfDNA isolation on the EZ1 Advanced XL instrument provides convenient, fully automated bead-based purification of up to 14 samples in parallel, from up to 8 ml of sample. cfDNA isolation on the EZ2 Connect instrument provides convenient, fully automated bead-based purification of up to 24 samples in parallel, from up to 8 ml of sample.

The eluted ccfDNA is ready for use in downstream reactions or storage at –30 to –15°C. Purified nucleic acids are free of proteins, nucleases, and other impurities.


Principle and procedure

The EZ1&2 ccfDNA Kit combines the speed and efficiency of anion exchange-based nucleic acid purification with the convenient automated handling of magnetic particles. The purification procedure comprises 3 steps: bind, wash, and elute. It is suited for simultaneous processing of multiple samples and provides pure nucleic acid in a short time for up to 14 samples on the EZ1 Advanced XL or up to 24 samples on the EZ2 Connect. Run times depending on the sample volume are as follows: 2 ml: 35 min, 4 ml: 46 min, 8 ml: 70 min. To allow processing of these volumes of plasma or serum, the kit contains Large Volume Tubes. To accommodate the Large Volume Tubes (7 ml), special Tip Racks have to be used. See details in section Automation and the Ordering Information (page 33).

To isolate circulating, cell-free nucleic acids from blood samples, plasma needs to be prepared from whole blood. The obtained plasma can be stored at -30°C to -15°C in the supplied Large Volume Tubes (7 ml) until the start of the automated ccfDNA isolation. For ccfDNA isolation, the plasma samples are loaded in the Large Volume Tubes onto the instrument.

Adsorption to anion exchange magnetic particles

Proteinase K and binding buffer are automatically added to the samples provided in the Large Volume Tubes to adjust binding conditions. Samples are thoroughly mixed with magnetic particles to allow adsorption of ccfDNA to the surface. Salt and pH conditions ensure that proteins and other contaminants, which can inhibit PCR and other downstream enzymatic reactions, are not bound to the magnetic particles.

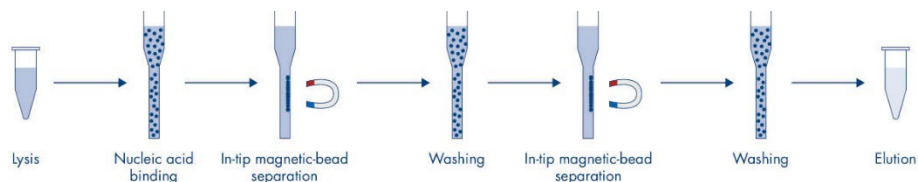

Washing of bound nucleic acids

While ccfDNA remains bound to the magnetic particles, contaminants are efficiently washed away during a series of 3 wash steps. Sample preparation waste is collected in tubes in positions 2, 3, and 4 of the tip rack, which holds the 7 ml large volume tubes to provide the sample and the tip holder. Remaining liquid in these positions is normal and not an indication of incomplete processing of the sample.

Elution of purified nucleic acids

Highly purified ccfDNA is eluted in low salt buffer (the eluate volume can range from 60–75 μl) and collected in 1.5 ml microcentrifuge tubes. The purified ccfDNA can be either used immediately in downstream applications or stored for future use.

We recommend storing purified circulating nucleic acids at 2–8°C if only for up to 24 hours or at –30 to –15°C if longer than 24 hours.

Automation

The EZ1&2 ccfDNA Kit can be used on the EZ1 Advanced XL instrument or the EZ2 Connect instrument. The provided protocol describes the workflow for both instruments. ccfDNA binding, washing, and elution steps are completely conducted on the EZ1 Advanced XL instrument or EZ2 Connect instrument.

To accommodate the Large Volume Tubes (7 ml) that are used with the EZ1&2 ccfDNA Kits, the EZ1 Advanced XL Tip Rack - Large Volume is required, which needs to be ordered separately (cat. no. 9027008) when using the EZ1 Advanced XL instrument. In addition, the EZ1 Advanced XL ccfDNA Card (cat. no. 9026964) providing the ccfDNA protocols is also required. When using the EZ2 Connect instrument, the EZ2 Connect Tip Rack - Large Volume is required, which needs to be ordered separately (cat. no. 9027011).

Starting material

The EZ1&2 ccfDNA Kit has been tested with serum and plasma samples.

Preparing plasma from whole blood

For recommendations for plasma preparation please refer to Appendix A: Recommendations for Plasma Separation and Storage”, page 27. This protocol includes a high *g*-force centrifugation step to remove cellular debris and reduce the amount of cellular or genomic DNA and RNA in the sample. The obtained plasma can be stored at –30°C to –15°C in the supplied Large Volume Tubes (7 ml) until the start of the automated ccfDNA isolation.

Sample volumes

The EZ1&2 ccfDNA Kit has been optimized for sample volumes 1–8 ml. Yield depends on the sample volume and the concentration of circulating nucleic acids in the sample (typically, 1–100 ng/ml in plasma). The EZ1&2 ccfDNA Kit has been designed to be flexible with serum or plasma input volumes and provides 3 different protocols to choose from depending on the input volume. The 2 ml protocol can be used to process sample volumes of up to 2 ml, the 4 ml protocol to process sample volumes of >2 ml to 4ml, and the 8 ml protocol to process samples >4 ml to 8 ml.

The provided large volume tubes can be used for storage of the generated plasma at –30 to –15°C until use on the instrument.

Sample numbers


The kit provides enough material to run 48 samples in total. Using the EZ1 Advanced XL instrument up to 14 samples can be run in parallel and up to 24 when using the EZ2 Connect.

Yield and size of nucleic acids

Yields of cell-free circulating nucleic acids isolated from biological samples are normally well below 1 µg and are therefore too low to determine with a spectrophotometer. Additionally, the small size of circulating DNA (main peak approximately 160–180 bp) renders many fluorometric methods unreliable for yield determination. Quantitative amplification methods, such as real-time PCR using a multi-copy target gene with known and stable copy number per genome (e.g., GAPDH or β-actin) is the most accurate method for determining the concentration of small amounts of cfDNA and recommended for determination of yields.

Size distribution of circulating nucleic acids purified using this procedure can be checked by analysis on an Agilent® Bioanalyzer®, TapeStation®, or a similar device.

EZ1&2 ccfDNA Workflow

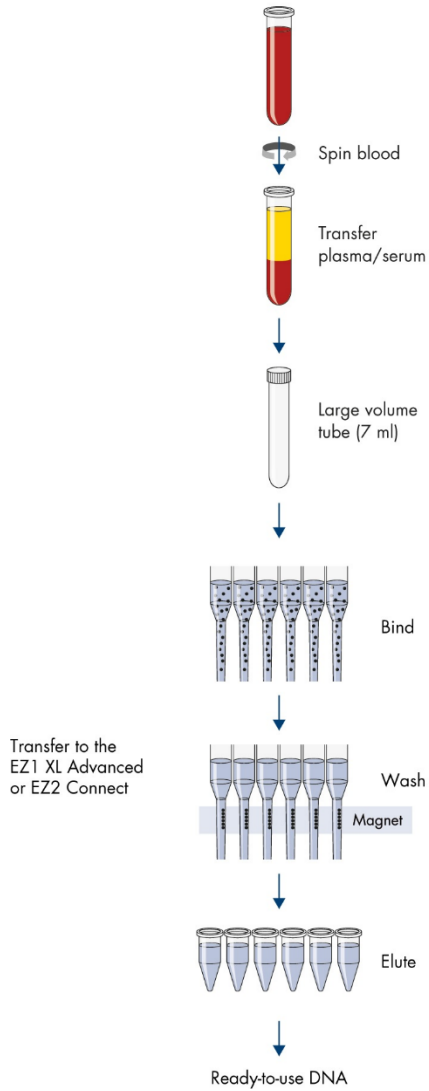

Equipment and Reagents to Be Supplied by User

When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, consult the appropriate safety data sheets (SDSs) available from the product supplier.

In addition to the EZ1&2 ccfDNA Kit, the following supplies are required.

- Vortex shaker for microcentrifuge tubes
- Pipettes (adjustable)
- Sterile pipette tips (pipette tips with aerosol barriers are recommended to help prevent cross-contamination)
- EZ1 Advanced XL:
 - EZ1 Advanced XL instrument (cat. no. 9001874)
 - EZ1 Advanced XL ccfDNA Card (cat. no. 9026964)
 - EZ1 Advanced XL Tip Rack - Large Volume (cat. no. 9027008)
- EZ2 Connect:
 - EZ2 Connect instrument (cat. no. 9003210)
 - EZ2 Connect Tip Rack - Large Volume (cat. no. 9027011)

Important Notes

- This handbook contains instructions on using the EZ1&2 ccfDNA Kit on the EZ1 Advanced XL and the EZ2 Connect instruments. The next section will describe the respective instruments and their specific requirements for this kit (e.g., use of a special rack, pages 16 and 19).
- The protocol on page 21 will include information on sample handling on the EZ1 Advanced XL (page 24) and the EZ2 Connect (page 23) instrument.
- UV decontamination helps to reduce possible pathogen contamination of the EZ1 Advanced XL worktable surfaces. The efficiency of inactivation has to be determined for each specific organism, and it depends, for example, on layer thickness and sample type. QIAGEN cannot guarantee complete inactivation of specific pathogens.

Reagent cartridges

Reagents for purification of nucleic acids from a single sample are contained in a single reagent cartridge (Figure 1). Each well of the cartridge contains a particular reagent, such as binding buffer or wash buffer.

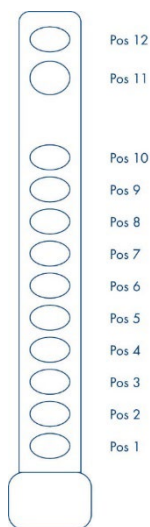
A**B**

Figure 1. Ease of worktable setup using reagent cartridges. A: A sealed, prefilled reagent cartridge (example). B: Loading reagent cartridges into the cartridge rack. The cartridge rack itself is labeled with an arrow to indicate the direction in which reagent cartridges must be loaded.

Working with the EZ2 Connect instrument

Worktable

The worktable of the EZ2 Connect instrument is where the user loads samples and the components of the EZ1&2 ccfDNA Kit (cat. no. 954854).

Details on worktable setup are provided in the protocol in this handbook and are also displayed on the touch screen when the user starts the worktable setup. The display also shows protocol status during the automated purification procedure.

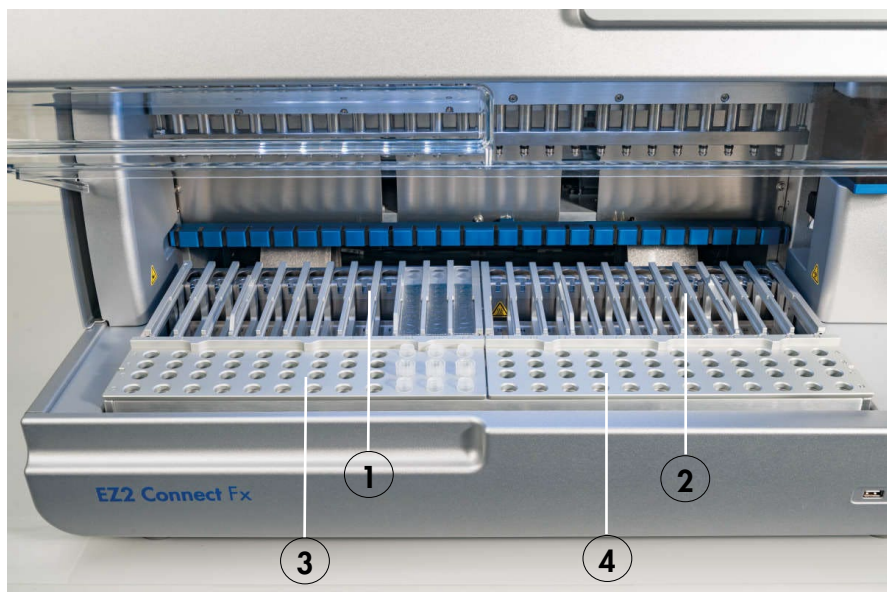


Figure 2. EZ2 Connect worktable.

- | | |
|--------------------------------------|---------------------------------------|
| 1. EZ2 Connect Cartridge Rack – left | 2. EZ2 Connect Cartridge Rack – right |
| 3. EZ2 Connect Tip Rack – left | 4. EZ2 Connect Tip Rack – right |

Operation of the EZ2 Connect

The EZ2 Connect provides various features to support the sample preparation workflow. These include functions for remote access via QIAsphere, data input via bar code reading, data storage and transfer, report generation, and guided instrument maintenance. For more information about these features, please refer to the *EZ2 Connect and EZ2 Connect Fx User Manual*.

EZ2 Connect Cartridge Rack

The EZ2 Connect Cartridge Rack consists of two separate parts. The left cartridge rack is used for cartridges in positions from 1 to 12. The right cartridge rack is used for cartridges in positions 13 to 24. On the worktable, the cartridge rack is located behind the tip rack (Figure 2). It holds up to 24 reagent cartridges.

EZ2 Connect Tip Rack



Figure 3. EZ2 Connect Tip Rack.

The EZ2 Connect Tip Rack consists of two separate parts. The left part of the tip rack is used for labware in positions from 1 to 12. The right part of the tip rack is used for labware in positions 13 to 24. The EZ2 Connect Tip Racks are located at the front of the worktable (Figure 2). Each consists of four rows and 12 columns. Individual positions in the EZ2 Connect Tip Rack are marked by engravings. During run setup, the user interface gives instructions to load specific positions of the EZ2 Connect Tip Rack with filter tips in tip holders, sample tubes, or elution tubes.

Important: To accommodate the Large Volume Tubes (7 ml) that are used with the EZ1&2 ccfDNA Kits, **the EZ2 Connect Tip Rack - Large Volume is required, which needs to be ordered separately (cat. no. 9027011)**. Remove the standard racks from the instrument and replace them with the EZ2 Connect Tip Rack - Large Volume.

Working with the EZ1 Advanced XL instrument

EZ1 Advanced XL Cards

Protocols for nucleic acid purification are stored on preprogrammed EZ1 Cards (integrated circuit cards). The user simply inserts an EZ1 Advanced XL Card into the EZ1 Advanced XL Instrument, and the instrument is then ready to run a protocol (Figure 4 and Figure 5). The availability of various protocols increases the flexibility of EZ1 instruments.



Figure 4. Ease of protocol setup using EZ1 Advanced XL Cards. Inserting an EZ1 Card, which contains a protocol, into the EZ1 Advanced XL instrument.



Figure 5. EZ1 Card completely inserted into EZ1 Card slot.

The instrument should only be switched on after an EZ1 Card is inserted. Make sure that the EZ1 Card is completely inserted (Figure 5). Otherwise, essential instrument data could be lost, leading to a memory error. Make sure to turn off the instrument before removing or replacing the EZ1 Cards; otherwise, a memory error can occur.

Worktable

The worktable of the EZ1 Advanced XL instrument is where the user loads samples and the components of the EZ1&2 ccfDNA Kit (Figure 6).

Details on worktable setup are provided in the protocol in this handbook and are also displayed in the vacuum fluorescent display (VFD) of the EZ1 Advanced XL control panel when the user starts the worktable setup. The display also shows protocol status during the automated purification procedure.

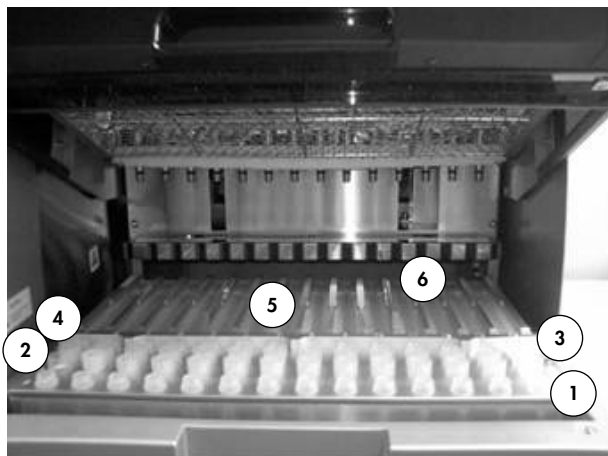


Figure 6. Worktable of an EZ1 Advanced XL Instrument.

- 1** First row: Elution tubes (ET) (1.5 ml) are loaded here.
- 2** Second row: Large volume tubes containing 2 or 4 ml of sample are loaded here. After the run is completed, the tubes contain liquid waste.
- 3** Third row: for 2 and 4 ml protocols, new (empty) large volume tubes are loaded here. For the 8 ml protocol, large volume tubes containing 4 ml of sample are loaded here. After the run is completed, the tubes contain liquid waste.
- 4** Fourth row: Disposable tip holders (DTH) containing disposable filter-tips (DFT) are loaded here.
- 5** Reagent cartridges (RCB) loaded into the cartridge rack.
- 6** The heating block positions contain magnetic bead suspension EZ and elution buffer EZE placed by the user before the start of the EZ1&2 ccfDNA Kit protocol.

EZ1 Advanced XL Tip Rack

Important: To accommodate the Large Volume Tubes (7 ml) that are used with the EZ1&2 ccfDNA Kits, **the EZ1 Advanced XL Tip Rack - Large Volume is required, which needs to be ordered separately (cat. no. 9027008) when using the EZ1 Advanced XL instrument.** Remove the standard racks from the instrument and replace them with the EZ2 Connect Tip Rack - Large Volume.

Data tracking with the EZ1 Advanced XL

The EZ1 Advanced XL enables complete tracking of a variety of data for increased process control and reliability. The EZ1&2 Kit lot number and expiration date are entered at the start of the protocol using the Q-Card bar code. A user ID and the Q-Card bar code can be entered manually using the keypad or by scanning bar codes using the handheld bar code reader. Sample and assay information can also be optionally entered at the start of the protocol. At the end of the protocol run, a report file is automatically generated. The EZ1 Advanced XL can store up to 10 report files, and the data can be transferred to a PC or directly printed using a printer.

To receive report files on a PC, the EZ1 Advanced Communicator software needs to be installed. The software receives the report file and stores it in a folder that you define. After the PC has received the report file, you can use and process the file with a LIMS (Laboratory Information Management System) or other programs. An example of a report file is shown in “Appendix B: Example of an EZ1 Advanced Report File”, page 30. In report files, the 14 pipetting channels of the EZ1 Advanced XL are named, from left to right, channels 1–14.

When scanning a user ID or Q-Card bar code with the bar code reader, a beep confirms data input. After the information is displayed for 2 seconds, it is automatically stored, and the next display message is shown. When scanning sample ID, assay kit ID, or notes, a beep confirms data input, the information is displayed, and a message prompts you to enter the next item of information. After scanning sample ID, assay kit ID, and notes, press **Ent** once to confirm that the information entered is correct. If, for example, a wrong bar code was scanned for one of the samples, press **Esc** and then rescan all sample bar codes according to the onscreen instructions. For user ID and notes, you can enter the numbers using the keypad, or you can easily generate your own bar codes to encode these numbers.

For details about data tracking and using EZ1 Advanced Communicator software, see the *EZ1 Advanced XL User Manual*, www.qiagen.com/HB-0176.

Protocol: Purification of ccfDNA from up to 8 ml Serum or Plasma Using the EZ1&2 ccfDNA Kit

Important points before starting

- If using the EZ1&2 ccfDNA Kit for the first time, read “Important Notes”, page 13.
- After receiving the kit, check the kit components for damage. If any kit components are damaged, contact QIAGEN Technical Services or your local distributor. In case of liquid spillage, refer to “Safety Information”, page 6. Do not use damaged kit components, because their use may lead to poor kit performance or contamination of the instrument.
- Make sure the EZ1 Advanced XL Tip Rack - Large Volume or EZ2 Connect Tip Rack - Large Volume is used instead of the regular tip rack.

Things to do before starting

- Plasma samples should be at room temperature before use. If plasma/serum was stored frozen, thaw the samples by incubating at 37°C in a water bath.
- If the sample volume is less than 2 ml, 4 ml, or 8 ml, adjust the volume with PBS to 2 ml, 4 ml, or 8 ml.
- It is very important to fill the large volume tubes with the volume specified (2 ml or 4 ml) and not more or less as this can lead to foaming or sample spillage resulting in poor ccfDNA quality and yield.
- EZ1 instruments should only be switched on after an EZ1 card is inserted. Make sure that the EZ1 Card is completely inserted; otherwise, essential instrument data could be lost, leading to a memory error. Make sure to turn off the instrument before removing or replacing the EZ1 Cards; otherwise, a memory error can occur.

Procedure

1. Load reagent cartridges into the cartridge rack.
2. Resuspend Magnetic Bead Suspension EZ by vortexing briefly for 1 min and then tap to ensure all the liquid is at the bottom of the tube. Remove the cap from the tube with the bead solution and place in position 11 of the EZ1&2 ccfDNA cartridge (see Figure 7).

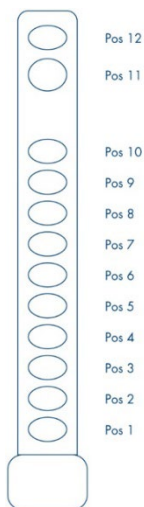


Figure 7. EZ1&2 ccfDNA cartridge.

3. Transfer 900 μ l elution buffer EZE into position 12 of the EZ1&2 ccfDNA cartridge (see Figure 7).
4. Remove caps of all tubes and prepare the large volume tip rack as follows (see Figure 8):
 - Position 4/A: Tip holder with Filter Tip (provided)
 - Position 3/B: new large volume tube (7 ml) (provided)
 - Position 2/C: new large volume tube (7 ml) (provided)
 - Position 1/D: new 1.5 ml elution tube (provided)

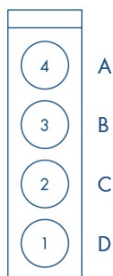


Figure 8. Tip Rack - Large Volume.

5. **2 ml/4 ml protocol:** Transfer exactly 2 ml/4 ml sample into the 7 ml tube in position 2. Leave the 7 ml tube in position 3 empty.

8 ml protocol: Split your sample and transfer exactly 4 ml into each of the 7 ml tubes in positions 2 and 3.

Procedure on the EZ2[®] Connect

1. Turn on the EZ2 Connect instrument.
2. Tap “DNA” on the Applications panel and select the “EZ1&2 ccfDNA Kit” and press **Next**. Follow onscreen instructions for selection of protocol, parameter definition, sample position selection, sample IDs, and worktable setup.
3. Open the instrument door. Load the cartridge rack into the instrument.
4. Place the EZ2 Connect Tip Rack - Large Volume into the instrument.
5. Close the instrument door; press **START** to initiate the EZ1&2 ccfDNA protocol.
6. The display will show “Protocol finished” when the run is completed. Select **Finish**.

Open the instrument hood. Remove the elution tubes containing the purified ccfDNA (in 75 µl) from position D of the tip rack. Discard the sample preparation waste (in tubes in positions 2 and 3 and in case of the 8 ml protocol in position 4) (see Figure 8).

Optional: Follow onscreen instructions for UV decontamination of worktable surfaces.

7. Perform regular maintenance after each run. Press **Finish** to return to the Home Screen.

Procedure on the EZ1 Advanced XL

1. For use on EZ1 Advanced XL instrument, insert the EZ1 Advanced XL ccfDNA Card completely into the card slot and turn on the instrument. Please make sure the instrument is turned off prior to inserting the card.
2. Press **START** to initiate the EZ1&2 ccfDNA protocol. Follow onscreen instructions for selection of protocol, data tracking, and worktable setup. Close the instrument door; press **START** to start the protocol.
3. Open instrument door. Load the cartridge rack into the instrument.
4. Place the EZ1 Advanced XL Tip Rack - Large Volume into the instrument.
5. The display shows "Protocol finished" when finished. Press **Esc**.
6. Open the instrument door. Remove elution tubes containing purified ccfDNA (in 75 μ l) from position 1 of the tip rack (see Figure 8). Discard the sample preparation waste (in tubes in positions 2 and 3 and in case of the 8 ml protocol in position 4) (see Figure 8). Press **Ent**. The report file is automatically transferred.
Optional: Follow onscreen instructions for UV decontamination of table surfaces.
7. Perform regular maintenance after each run. Press **Esc** to return to the main menu.
Note: Regular maintenance consists of cleaning the piercing unit and the worktable surfaces.
Important: The piercing unit is sharp. The use of double gloves is recommended.
8. To run another protocol, press **Start** and follow the protocol from step 1. Otherwise, close the instrument door and switch off the instrument.

Troubleshooting Guide

This troubleshooting guide may be helpful in solving any problems that may arise. For more information, see also the Frequently Asked Questions page at our Technical Support Center: www.qiagen.com/FAQ/FAQList.aspx. The scientists in QIAGEN Technical Services are always happy to answer any questions you may have about either the information and protocols in this handbook or sample and assay technologies (for contact information, see back cover or visit www.qiagen.com).

Comments and suggestions

Little or no nucleic acids in the eluate

- | | |
|---|---|
| a) Samples frozen and thawed more than once | Repeated freezing and thawing should be avoided because this may lead to DNA degradation. Always use fresh samples or samples that have been thawed only once. |
| b) Low concentration of target DNA in the samples | Samples were left standing at room temperature for too long. Repeat the purification procedure with new samples. |
| c) Magnetic particles not completely resuspended | Make sure to resuspend the magnetic particles thoroughly and no beads remain in the caps when they are removed, before loading the tubes into position 11 of the reagent cartridges. |
| d) Foaming during purification | Make sure that the large volume tubes are filled with 2 ml (2ml protocol) or 4 ml (4 ml or 8 ml protocol for both tubes) of sample and not less as this can cause foaming, which may drastically reduce ccfDNA yield and purity. |
| e) Wrong measurement method used | Due to low yields and the small size of ccfDNA fragments, not all measurement methods are reliable, see "Yield and size of nucleic acids", page 10. |
| f) Elution buffer EZE expired | Exchange with surrounding air may lead to reduced stability of elution buffer EZE leading to reduced efficiency of ccfDNA extraction. Close partially used bottles carefully after use, and do not use remaining elution buffer EZE more than 4 weeks after opening the bottle. The volume of elution buffer EZE provided is sufficient to allow discarding the last 2–3 ml of each bottle. |

Comments and suggestions

Co-purification of cellular DNA

- | | |
|--|--|
| a) Extended time between blood draw and plasma preparation | Blood cells may disintegrate and release genomic DNA into the plasma, diluting the target nucleic acid. If blood cannot be processed to plasma immediately, use of dedicated collection tubes with ccfDNA stabilization, e.g., PAXgene Blood ccfDNA Tubes, is recommended. |
|--|--|

DNA does not perform well in downstream enzymatic reactions

- | | |
|--|---|
| a) Little or no DNA in the eluate | See "Little or no nucleic acids in the eluate" above for possible reasons. Increase the amount of eluate added to the reaction if possible. |
| b) New <i>Taq</i> DNA polymerase or PCR chemistry used | If enzymes are changed, it may be necessary to readjust the amount of eluate used for PCR. |

General handling

- | | |
|--|---|
| a) Error message in instrument display | Refer to the user manual supplied with your EZ1 Advanced XL or EZ2 Connect instrument. |
| b) Report file not printed | Check whether the printer is connected to the EZ1 Advanced XL via the "PC/Printer" serial port. Check whether the serial port is set for use with a printer. |
| c) Wrong Q-Card ID entered | If the wrong ID was entered instead of the Q-Card ID, the EZ1 Advanced XL will not accept the ID and will prompt for the Q-Card ID until the correct ID is entered. Press Stop twice to go to the main menu. |

References

1. Sambrook, J., and Russell, D.W. (2001). *Molecular Cloning: A Laboratory Manual*. 3rd ed. Cold Spring Harbor, NY: Cold Spring Harbor Laboratory Press).

Appendix A: Recommendations for Plasma Separation and Storage

To isolate circulating, cell-free nucleic acids from blood samples, we recommend following this protocol, which includes a high g -force centrifugation step to remove cellular debris and reduce the amount of cellular or genomic DNA and RNA in the sample. Human serum or plasma samples can be generated using blood collection tubes such as the PAXgene Blood ccfDNA Tube; please refer to the manufacturer's recommendation for plasma separation procedure.

Procedure

1. Place whole blood in a collection tube into a centrifuge with a swing-out rotor and appropriate buckets.
2. Centrifuge the blood samples for 10 min at $1,900 \times g$ with temperature set to 4°C .
3. Carefully aspirate plasma supernatant without disturbing the buffy coat layer. Approximately 4–5 ml plasma can be obtained from one 10 ml primary blood tube.

Note: Plasma can be used for circulating nucleic acid extraction at this stage. However, the following high-speed centrifugation will remove additional cellular debris and contamination of the circulating nucleic acids by genomic DNA and RNA derived from damaged blood cells.

4. Transfer aspirated plasma into new 15 ml centrifuge tubes with conical bottoms.
5. Centrifuge the plasma samples for 15 min at $3,000 \times g$ in a fixed-angle rotor with temperature set to 4°C .

This will remove additional cellular nucleic acids attached to cell debris.

Note: For certain blood collection tubes, specific centrifugation conditions may be recommended by the manufacturer, e.g., for the PAXgene Blood ccfDNA tubes, for the second centrifugation, 15 min at $1,900 \times g$ is recommended.

The second centrifugation may also be performed at 10,000 or up to 16,000 $\times g$. This will remove some of the larger extracellular vesicles (EVs) and cell fragments, which may also contain cell-free DNA.

6. Using a pipette, carefully transfer the supernatant into a new tube without disturbing the pellet.
7. If plasma will be used for nucleic acid extraction on the same day, store at 2–8°C until further processing. For longer storage, keep plasma frozen at –90 to –65°C. Before using the plasma for circulating nucleic acid extraction, thaw plasma tubes at room temperature or at 37°C in a water bath.
8. In case of cryoprecipitates, follow these 2 steps:
 - 8a. Centrifuge the plasma sample using the same conditions as in step 5.
 - 8b. Transfer the supernatant into a new tube, and then begin with the nucleic acid extraction protocol.

Appendix B: Example of an EZ1 Advanced Report File

This appendix shows a typical report file generated on the EZ1 Advanced. The values for each parameter will differ from the report file generated on your EZ1 Advanced.

Please note that **User ID** is allowed a maximum of 9 characters, and that **Assay kit ID** and **Note** are allowed a maximum of 14 characters.

The EZ1 Advanced XL generates a similar report file containing instrument and protocol information relevant to the EZ1 Advanced XL and information for channels 1–14.

REPORT - FILE EZ1 Advanced:

```
-----  
Serial no. EZ1 Advanced: _____ SN 0001  
User ID: _____ 9876543210  
Firmware version: _____ V 1.0.0  
Installation date of instrument: _____ Jan 05, 2008  
Weekly maintenance done on: _____ Jun 15, 2008  
Yearly maintenance done on: _____ Jan 10, 2008  
Date of last UV-run: _____ Mar13, 2008  
Start of last UV-run: _____ 16:06  
End of last UV-run: _____ 16:26  
Status UV-run: _____ o.k.  
Protocol name: _____ Virus 2.0  
-----  
Date of run: _____ Jul 25, 2008
```

Start of run: _____ 12:57

End of run: _____ 13:50

Status run: _____ o.k.

Error code: _____

Sample input volume [ul]: _____ 400

Elution volume [ul]: _____ 150

Channel A:

Sample ID: _____ 123456789

Reagent kit number: _____ 9801401

Reagent lot number: _____ 1181234567

Reagent expiry date: _____ 1210

Assay kit ID: _____ 848373922

Note: _____ 2000

EZ1 Virus Handbook 10/2010 41

Channel B:

Sample ID: _____ 234567890

Reagent kit number: _____ 9801401

Reagent lot number: _____ 1181234567

Reagent expiry date: _____ 1210

Assay kit ID: _____ 836266738

Note: _____

Channel C:

Sample ID: _____ 345678901

Reagent kit number: _____ 9801401

Reagent lot number: _____ 1181234567

Reagent expiry date: _____ 1210

Assay kit ID: _____ 883727832

Note: _____ 1000

Channel D:

Sample ID: _____ 456789012

Reagent kit number: _____ 9801401

Reagent lot number: _____ 1181234567

Reagent expiry date: _____ 1210

Assay kit ID: _____ 763684837

Note: _____

Channel E:

Sample ID: _____ 567890123

Reagent kit number: _____ 9801401

Reagent lot number: _____ 1181234567

Reagent expiry date: _____ 1210

Assay kit ID: _____ 4387728002

Note: _____

Channel F:

Sample ID: _____ 678901234

Reagent kit number: _____ 9801401

Reagent lot number: _____ 1181234567

Reagent expiry date: _____ 1210

Assay kit ID: _____ 509389403

Note: _____ 50

Ordering Information

Product	Contents	Cat. no.
EZ1&2 ccfDNA Kit (48)	48 reagent cartridges (EZ1&2 ccfDNA), 96 Large Volume Tubes 7 ml, Elution Buffer EZE, Magnetic Bead Suspension EZ, Buffers, 50 Elution Tubes, 50 Disposable Tip Holders, 50 Disposable Filter-Tips, 100 Elution Tubes (1.5 ml)	954854
EZ1 Advanced XL Tip Rack - Large Volume	Metal rack that accommodates tip holders containing filter-tips, elution tubes, and large volume tubes on the worktable.	9027008
EZ1 Advanced XL, System	Robotic workstation for automated purification of nucleic acids from up to 14 samples using EZ1 Kits: includes installation, training, 1-year warranty on parts and labor	9001874
EZ1 Advanced XL ccfDNA Card	Preprogrammed card for EZ1 ccfDNA protocol on the EZ1 Advanced XL instrument	9026964
EZ2 Connect instrument	Benchtop instrument for automated isolation of nucleic acids from up to 24 samples in parallel, using sealed prefilled cartridges; includes 1-year warranty on parts and labor.	9003210
EZ2 Connect Tip Rack - Large Volume	Metal rack that accommodates tip holders containing filter-tips, elution tubes, and large volume tubes on the worktable.	9027011

Product	Contents	Cat. no.
Related products		
PAXgene Blood ccfDNA Tubes (100), RUO*	100 blood collection tubes (10 ml)	768115
Accessories		
Filter Tips and Holders, EZ1 (50)	50 Disposable Filter-Tips, 50 Disposable Tip Holders; additional tips and holders for use with EZ1 Kits	994900

* Not available in all countries

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual. QIAGEN kit handbooks and user manuals are available at www.qiagen.com or can be requested from QIAGEN Technical Services or your local distributor.

Document Revision History

Date	Changes
01/2022	Initial version.

Notes

Limited License Agreement for EZ1&2 ccfDNA Kit

Use of this product signifies the agreement of any purchaser or user of the product to the following terms:

1. The product may be used solely in accordance with the protocols provided with the product and this handbook and for use with components contained in the kit only. QIAGEN grants no license under any of its intellectual property to use or incorporate the enclosed components of this kit with any components not included within this kit except as described in the protocols provided with the product, this handbook, and additional protocols available at www.qiagen.com. Some of these additional protocols have been provided by QIAGEN users for QIAGEN users. These protocols have not been thoroughly tested or optimized by QIAGEN. QIAGEN neither guarantees them nor warrants that they do not infringe the rights of third-parties.
2. Other than expressly stated licenses, QIAGEN makes no warranty that this kit and/or its use(s) do not infringe the rights of third-parties.
3. This kit and its components are licensed for one-time use and may not be reused, refurbished, or resold.
4. QIAGEN specifically disclaims any other licenses, expressed or implied other than those expressly stated.
5. The purchaser and user of the kit agree not to take or permit anyone else to take any steps that could lead to or facilitate any acts prohibited above. QIAGEN may enforce the prohibitions of this Limited License Agreement in any Court, and shall recover all its investigative and Court costs, including attorney fees, in any action to enforce this Limited License Agreement or any of its intellectual property rights relating to the kit and/or its components.

For updated license terms, see www.qiagen.com.

Trademarks: QIAGEN[®], Sample to Insight[®], EZ1[®] (QIAGEN Group); Agilent[®], Bioanalyzer[®], TapeStation[®] (Agilent Technologies, Inc.); PAXgene[®] (PreAnalytiX GmbH). Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

01/2022 HB-2915-002 © 2022 QIAGEN, all rights reserved.

Ordering www.qiagen.com/shop | Technical Support support.qiagen.com | Website www.qiagen.com