

Data Sheet

dPCR CGT Assay

dPCR CGT assays	bGH poly A (FAM, HEX, Cy5)
Assay catalog numbers	250233, 250234, 250235
dPCR instrument	QIAcuity Digital PCR System
dPCR master mix	QIAcuity Probe PCR Kit
Probe fluorophores and quencher	FAM, HEX: ZEN/Iowa Black Cy5: TAO/Iowa Black
Amplicon length	70 bp
Restriction enzyme compatibility	<i>MspI</i> , <i>HpaII</i> , <i>SmaI</i> , <i>HaeIII</i>

Description


The bGH poly A assay is a wet-lab validated, 20x assay targeting the bovine growth hormone polyadenylation signal sequence.

Synthetic DNA sequence for positive controls

ATATCTCGACTGTGCCTTCTAGTTGCCAGCCATCTGTTGTTTGCCCCCTCCCCCGTGCCTTCC
TTGACCCTGGAAGGTGCCACTCCCACTGTCCTTTCCTAATAAAATGAGGA

Sample data



1D scatter plot: Single-well data of ITR assays (Green, Yellow and Crimson channel) on lysed adeno-associated virus (AAV) sample.

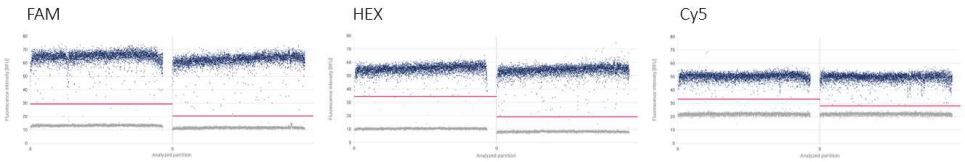

2D scatter plot: Single-well data of (A) bGH poly A assay (Crimson channel, Y-axis) and CMV promoter assay (Yellow channel, X-axis) or (B) bGH poly A assay (Green channel, X-axis) and ITR assay (Yellow channel, Y-axis) or (C) bGH poly A assay (Yellow channel, X-axis) and CMV promoter assay (Crimson channel, Y-axis) with AAV2 DNA input.

A

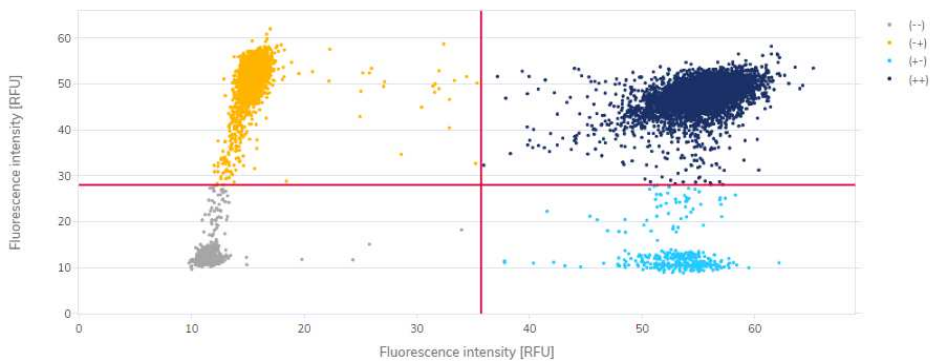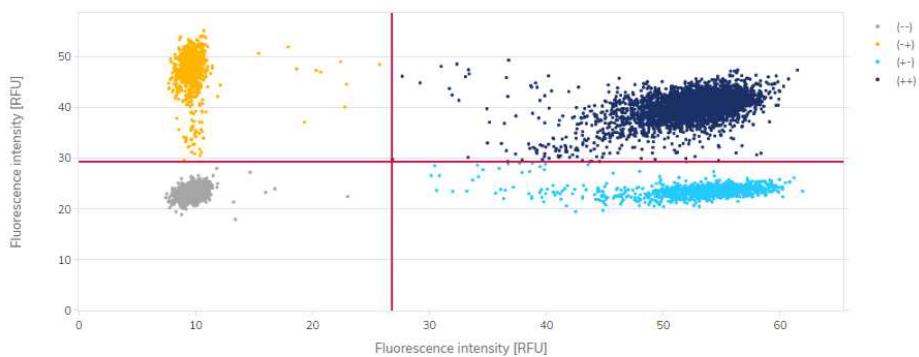

B**C**

Document Revision History

Date	Changes
06/2022	Initial release

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual.

Trademarks: QIAGEN®, Sample to Insight®, QIAcuity® (QIAGEN Group). Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

06/2022 HB-3080-001 © 2022 QIAGEN, all rights reserved.